

# College Catalog

2025

# Welcome!

Dear Educators and Partners,

What a warm welcome I have received this year! I'm so thrilled to join the mission that Dave and Christine created at Vernier and to work alongside educators who are at the heart of everything we do. My decision to join Vernier as CEO was deeply personal—shaped by my father's career as a teacher and administrator—and grounded in a shared commitment to advancing science education through hands-on learning and tools that make STEM accessible and engaging for all students.

At Vernier, we know that STEM education is transformative. Our solutions are not just designed for data collection but to empower students to become active participants in their learning. When students can observe real-world phenomena, ask questions, conduct investigations, and analyze data, they move beyond what's possible from only reading or watching to gain a deeper understanding of science. This experiential approach helps students build critical thinking skills and fosters a resilient mindset that will serve them in any path they choose, from continuing in higher education to entering the workforce.

Vernier data-collection tools transform a classroom from a place where science is talked about to a place where it is actively experienced. Through real-time data collection and analysis, students begin to form critical connections and see firsthand how science shapes the world around them. This is where the inspiration happens—where science becomes personal, exciting, and meaningful. Our goal is for students to see themselves as real scientists and problem solvers, gaining confidence to navigate our rapidly changing world with knowledge, curiosity, and skill.

For me, the chance to support educators and create these moments of inspiration for students is a privilege and a passion. I am excited to continue the Vernier legacy of excellence and to work alongside all of you as we inspire the next generation of thinkers, explorers, and problem solvers. Together, let's make science accessible, engaging, and empowering for every student.

Warm regards,



A handwritten signature in black ink, appearing to read 'J Hedrick', positioned to the right of the circular portrait.

**Jill Hedrick**  
CEO

## Why Vernier? Endless Possibilities.

Our durable hardware and quality software are designed and priced for hands-on student use whether learning remotely or in the laboratory. We have ready-to-go experiments and resources in a wide variety of subjects.

BIOLOGY • CHEMISTRY • PHYSICS • ENGINEERING  
AGRICULTURAL SCIENCE • ENVIRONMENTAL SCIENCE • PHYSIOLOGY

Our sensors and data-collection technology are so versatile that you can use them in nearly any science or engineering course.

## What's New

### Accessories for an Eco-Friendly Lab

Three new products are now a part of the Vernier college chemistry lineup, offering reliable, easy-to-use solutions that enhance the lab experience for students and educators. These accessories can be used alongside Vernier Go Direct® sensors to provide impactful learning experiences that are both efficient and environmentally conscious.

Learn more on page 16.



ELECTROCHEMISTRY HALF-CELL PLATE  
ECHEM-PLT \$59

[vernier.com/echem-plt](http://vernier.com/echem-plt)



ELECTROCHEMISTRY METALS KIT  
ECHEM-MTLS ⬆ \$20

[vernier.com/echem-mtls](http://vernier.com/echem-mtls)



PRECISION VOLUME DISPENSER  
DC-DISP \$85

[vernier.com/dc-disp](http://vernier.com/dc-disp)

## Now Shipping

### Go Direct Salinity

Go Direct® Salinity precisely measures the total dissolved salt content of ocean or brackish water.

GDX-SAL \$180

[vernier.com/gdx-sal](http://vernier.com/gdx-sal)



### Go Direct Soil Moisture

This sensor uses capacitance to measure the volumetric water content of soil.

GDX-SM \$149

[vernier.com/gdx-sm](http://vernier.com/gdx-sm)



Watch the video!

### Go Direct Force Plate

Measure the forces developed during stepping, jumping, and other human-scale actions. Observe change in apparent weight as you ride an elevator or measure reaction forces as you lean against a wall.

GDX-FP \$399

[vernier.com/gdx-fp](http://vernier.com/gdx-fp)



Watch the video!

### Cart Fan

Add this modular fan to the Dynamics Cart and Track System and Go Direct Sensor Cart to provide a constant force for kinematics and dynamics experiments.

DTS-CFAN \$89

[vernier.com/dts-cfan](http://vernier.com/dts-cfan)



Watch the video!



Students don't learn effectively by just sitting and listening—they need to be consistently engaged in what they are learning about. As educators, we are always looking for ways to optimize the learning experience for our students, especially in today's increasingly remote environment.

*David Craig*  
*Oregon State University*

## Join these institutions, and hundreds of others, already using Vernier technology:

- |                                       |                                      |                                     |
|---------------------------------------|--------------------------------------|-------------------------------------|
| Arizona State University              | McGill University                    | University of Cambridge             |
| Baltimore City Community College      | Miami University                     | University of Chicago               |
| Benedictine University                | Michigan Technological University    | University of Hong Kong             |
| California State University—Fullerton | Mississippi State University         | University of Kansas                |
| Cameron University                    | National University of Colombia      | University of Minnesota—Minneapolis |
| Canisius College                      | National University of Singapore     | University of Nebraska—Lincoln      |
| Charles University                    | Oregon State University              | University of Pennsylvania          |
| Colorado School of Mines              | Princeton University                 | University of Puerto Rico           |
| Cornell University                    | Queensborough Community College      | University of Sydney                |
| Cuyahoga Community College            | Quinnipiac University                | University of Tennessee—Chattanooga |
| Delft University of Technology        | San José State University            | University of Toronto               |
| Dickinson College                     | Stanford University                  | University of Washington—Seattle    |
| ETH Zurich                            | Stephen F. Austin State University   | University of Wisconsin—Madison     |
| Georgia Tech                          | Sungkyunkwan University              | Vincennes University                |
| Harvard University                    | TEC Monterrey                        | Virginia Commonwealth University    |
| Haskell Indian Nations University     | Texas A&M                            | Wake Technical Community College    |
| Immaculata University                 | The Ohio State University            | West Virginia Wesleyan College      |
| Lehigh University                     | University of Arizona                | Yale University                     |
| Lund University                       | University of British Columbia       |                                     |
| Massachusetts Institute of Technology | University of California—Los Angeles |                                     |

## Partnerships

### LabArchives

We have partnered with LabArchives™ to bring high-quality biology content to instructors through the Lab Builder library. Because all content is structured and standardized, instructors can arrange, customize, and add content to their courses with ease. [Learn more at vernier.com/lab-archives](https://www.vernier.com/lab-archives)

### ADInstruments

We have partnered with ADInstruments to provide state-of-the-art, customizable biology content using ADInstruments' platform, Lt. This platform supports real-time data collection and analysis with many Vernier sensors and includes instructional videos, quizzes, and classroom management tools. [Learn more at vernier.com/adinstruments](https://www.vernier.com/adinstruments)

# Biology

## Why Vernier?

Vernier biology solutions help students form a deep understanding of key scientific concepts. Whether you are introducing your students to enzymes or exploring primary productivity, our probeware and ready-to-go experiments are the right fit for your laboratory.

### Quality

Durable hardware for lab and field use

### Affordable

Designed for education and education budgets

### Versatile

Supports a variety of devices and experiments



“

Your great products and superb support of them have been a major part of my labs and are very much appreciated.

*David Willey  
University of Pittsburgh*

[vernier.com/biology](https://www.vernier.com/biology)

# A Guide to Vernier Data Collection

We're here to support you as an educator as you incorporate data-collection technology into your instruction. See how our products provide you with affordable laboratory solutions designed for student success.

Our Guarantee: Most of our products are protected by a 5-year limited warranty. And after 5 years? We'll make every attempt to repair your equipment.

## What You Need to Get Started with Go Direct Sensors

### Go Direct Sensor

These versatile sensors connect to your device via Bluetooth® wireless technology or USB.

### Device

Go Direct® sensors connect to a wide variety of commonly used devices, including Chromebooks, computers, tablets, smartphones, and LabQuest® 3.

### Software

Vernier Graphical Analysis® Pro  
Vernier Spectral Analysis®  
Vernier Instrumental Analysis®

### Lab Book

Our popular, award-winning lab books provide hundreds of well-tested, customizable experiments. Our lab books come with a generous site license—purchase once and share files across your department.

## What You Need to Get Started with LabQuest 3

### Sensor

**Go Direct Sensor**  
These versatile sensors connect to LabQuest 3 via Bluetooth wireless technology or USB.  
*Go Direct spectrometers connect only via USB.*

**LabQuest Sensor**  
LabQuest sensors connect directly to LabQuest 3 sensor ports (BTA/BTD).

### LabQuest 3

LabQuest 3 serves as a standalone data-collection platform that works with all Vernier sensors.

### Software

LabQuest App

### Lab Book

Our popular, award-winning lab books provide hundreds of well-tested, customizable experiments. Our lab books come with a generous site license—purchase once and share files across your department.

## Software

### Vernier Graphical Analysis Pro

Our award-winning app enables advanced data analysis, provides opportunities to reinforce and extend learning, and includes features that support hybrid teaching models. We're always working to add and improve features—check out the “What's New” tab in the app.

**1-Year Site License (unlimited seats)** GAP-1YR \$149

**3-Year Site License (unlimited seats)** GAP-3YR \$299

#### Free Trial for Educators

Try out Graphical Analysis Pro for free for 30 days. Access the sample experiments and enhanced analysis tools to use with your students. Get a free trial and learn about site license options at [vernier.com/graphical-analysis](https://www.vernier.com/graphical-analysis)



### Vernier Spectral Analysis

Spectral Analysis supports our family of spectrometers on computers, Chromebooks, and compatible mobile devices. Use it to generate full spectra, create standard curves, and conduct kinetics experiments.

[vernier.com/spectral-analysis](https://www.vernier.com/spectral-analysis)

### Vernier Instrumental Analysis

Vernier Instrumental Analysis is used for more advanced instrumentation such as Go Direct Mini GC,™ Go Direct Polarimeter, and Go Direct Cyclic Voltammetry System.

[vernier.com/instrumental-analysis](https://www.vernier.com/instrumental-analysis)

Get everything you need with our new 4-in-1 software bundle!  
**3-Year Site License** DCA-3YR \$499 | Learn more at [vernier.com/dca-3yr](https://www.vernier.com/dca-3yr)

## LabQuest 3



**LabQuest 3 is a powerful, advanced, easy-to-navigate, and versatile data-logging solution for STEM students.**

LabQuest 3 is a standalone data-collection platform that students can use to collect, analyze, and interact with data efficiently. With its touch-screen abilities, students can navigate the platform with ease, and because of its wireless capabilities, students can collect data anywhere.

LABQ3 \$399

[vernier.com/labq3](https://www.vernier.com/labq3)

### LabQuest App

LabQuest 3 has built-in software that gives your students real-time graphing capabilities in a handheld device. It's powerful, yet beautifully simple.



# Free 30-Day Preview

Not sure if our technology is right for your laboratory? Our risk-free preview program makes it easier to decide. Most Vernier products are available for a 30-day preview (or longer, if requested) to educational institutions.

## General Biology

### Go Direct CO<sub>2</sub> Gas

This sensor measures gaseous carbon dioxide concentration levels, air temperature, and relative humidity. With built-in temperature compensation and humidity protection, this sensor is ideal for measuring fermentation, respiration, and photosynthesis rates.

GDX-CO2 \$229

[vernier.com/gdx-co2](http://vernier.com/gdx-co2)



### Go Direct SpectroVis Plus

Use Go Direct SpectroVis® Plus to collect a full-wavelength spectrum (absorbance, percent transmittance, fluorescence, or intensity), study absorbance vs. concentration (create standard curves), or monitor enzyme activity (enzyme kinetics).

GDX-SVISPL \$479

[vernier.com/gdx-svispl](http://vernier.com/gdx-svispl)



**Vernier Spectral Analysis**  
FREE DOWNLOAD  
Learn more on page 7.

### Go Direct Tris-Compatible Flat pH

Use this sensor to measure the pH of solutions. It features a sealed, gel-filled, double-junction electrode, making it compatible with Tris buffers and solutions containing proteins or sulfides.

GDX-FPH \$139

[vernier.com/gdx-fph](http://vernier.com/gdx-fph)



### Go Direct Optical Dissolved Oxygen

Use this sensor to measure dissolved oxygen, water temperature, and atmospheric pressure. It's ideal for experiments in biology, ecology, and environmental science.

GDX-ODO \$385

[vernier.com/gdx-odo](http://vernier.com/gdx-odo)




### Biology with Vernier

This book includes 31 experiments for fundamental concepts in biology. The instructor information section included for each experiment contains reagent preparation information, sample data, and tips for successful completion.

#### Topics

- Cell respiration
- Membrane diffusion
- Osmosis
- Photosynthesis and transpiration
- Human physiology

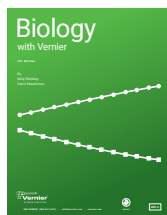

**Printed book + download**

BWV \$59

**Download only**

BWV-E \$45


[vernier.com/bwv](http://vernier.com/bwv)

### Investigating Biology through Inquiry

This book includes 22 investigations for many fundamental concepts in biology. Each investigation includes a preliminary activity, instructor information, sample researchable questions, and sample data.

#### Topics

- Cell and molecular biology
- Organismal biology
- Ecology
- Evolution

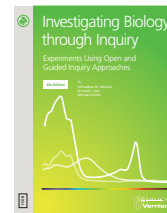

**Printed book + download**

BIO-I \$59

**Download only**

BIO-I-E \$45

[vernier.com/bio-i](http://vernier.com/bio-i)

### Biology Go Direct Standard Package

This package includes 11 sensors that work with Vernier Graphical Analysis Pro and LabQuest 3.

- Go Direct Temperature Probe
- Go Wireless® Heart Rate
- Go Direct Gas Pressure Sensor
- Go Direct CO<sub>2</sub> Gas Sensor
- Go Direct O<sub>2</sub> Gas Sensor
- Go Direct Conductivity Probe

- Go Direct pH Sensor
- Go Direct Colorimeter
- Go Direct Optical Dissolved Oxygen Probe
- Go Direct EKG Sensor
- BioChamber 250 ⚡
- BioChamber 2000 ⚡
- Go Direct Respiration Belt

GDP-BIO-ODX ⚡ \$1,820

[Learn more at vernier.com/gdp-bio-odx](http://vernier.com/gdp-bio-odx)



Buy 8 or more packages at \$1,784 and save \$288.

⚡ For California Proposition 65 warning, see page 36.

# Human Physiology

## Go Direct EKG

Use Go Direct® EKG to record electrical activity of the heart or skeletal muscles.

GDX-EKG \$175 [vernier.com/gdx-ekg](http://vernier.com/gdx-ekg)



## Go Direct Hand Dynamometer

Measure grip and pinch strength, and perform muscle fatigue studies.

GDX-HD \$135 [vernier.com/gdx-hd](http://vernier.com/gdx-hd)



## Go Direct Respiration Belt

Use this sensor to measure human respiration rate and study breathing patterns.

GDX-RB \$125

[vernier.com/gdx-rb](http://vernier.com/gdx-rb)



## Go Direct Spirometer

This multi-channel sensor can be used to measure tidal volume, vital capacity, flow rate, air pressure, and respiration rate.

GDX-SPR \$229

[vernier.com/gdx-spr](http://vernier.com/gdx-spr)



## Go Direct O<sub>2</sub> Gas

Use this sensor to measure gaseous oxygen concentration levels and air temperature.

GDX-O2 \$219

[vernier.com/gdx-o2](http://vernier.com/gdx-o2)



## Go Direct Blood Pressure

This affordable, non-invasive sensor is designed to easily measure human blood pressure.

GDX-BP \$125 [vernier.com/gdx-bp](http://vernier.com/gdx-bp)



## Human Physiology Experiments: Volume 1

This book contains 14 experiments that encourage students to investigate the physiology of the cardiac, muscular, respiratory, vascular, and nervous systems using Go Direct sensors.

[vernier.com/hsb-hp](http://vernier.com/hsb-hp)

*This lab book provides instructions for data collection with Vernier Graphical Analysis® only.*



**Printed book + download**

HSB-HP \$45

**Download only**

HSB-HP-E \$35

## Human Physiology Experiments: Volume 2

An expansion of our *Human Physiology Experiments: Volume 1* lab book, the setup for these experiments is minimal—students are collecting data within minutes.

[vernier.com/alb-hp2](http://vernier.com/alb-hp2)

*This lab book provides instructions for data collection with Vernier Graphical Analysis only.*



**Printed book + download**

ALB-HP2 \$45

**Download only**

ALB-HP2-E \$35

## Human Physiology Go Direct Standard Package

This package includes 9 sensors that work with Vernier Graphical Analysis Pro and LabQuest® 3. Two useful accessories are also included.

### SENSORS

- Go Direct EKG
- Go Direct Force and Acceleration
- Go Direct Surface Temperature
- Go Direct Hand Dynamometer
- Go Direct Respiration Belt
- Go Direct O<sub>2</sub> Gas
- Go Direct Blood Pressure
- Go Direct Spirometer
- Go Wireless® Heart Rate

### ACCESSORIES

- Reflex Hammer Accessory Kit
- BioChamber 250 ♣

GDP-HP-DX ♣ \$1,363

Buy 8 or more packages at \$1,336 and save \$216.

**Starter package also available**

Learn more at [vernier.com/gdp-hp-dx](http://vernier.com/gdp-hp-dx)



♣ For California Proposition 65 warning, see page 36.



## Spectrometers

### Go Direct SpectroVis Plus

Use Go Direct SpectroVis® Plus to collect a full-wavelength spectrum (absorbance, percent transmittance, fluorescence, or intensity\*), study absorbance vs. concentration (create standard curves), or monitor enzymatic activity (enzyme kinetics).

GDX-SVISPL \$479

[vernier.com/gdx-svispl](http://vernier.com/gdx-svispl)



\* Fiber optic cable purchased separately

### Go Direct UV-VIS Spectrophotometer

The Go Direct UV-VIS Spectrophotometer connects to your device via Bluetooth® wireless technology or USB to generate full spectra, Beer's law data, and kinetic traces of ultraviolet and visible-absorbing samples such as aspirin, DNA, proteins, and NADH.

GDX-SPEC-UV \$2,999

[vernier.com/gdx-spec-uv](http://vernier.com/gdx-spec-uv)




### Go Direct Fluorescence/UV-VIS Spectrophotometer

This spectrophotometer measures the fluorescence and absorbance spectra of samples such as quinine sulfate, fluorescein, chlorophyll, and fluorescent proteins, all while connecting to your device via Bluetooth wireless technology or USB.

GDX-SPEC-FUV \$3,499

[vernier.com/gdx-spec-fuv](http://vernier.com/gdx-spec-fuv)

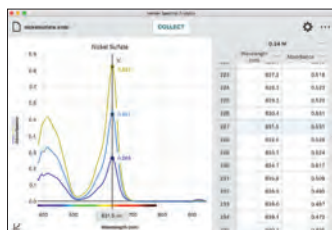

## Vernier Spectral Analysis

Our free Vernier Spectral Analysis® app makes it easy to incorporate spectroscopy into your general biology and biotechnology experiments. Using the app, students can collect a full spectrum and explore topics such as plant pigments, enzyme kinetics, and Beer's law (standard curves).

FREE DOWNLOAD [vernier.com/spectral-analysis](http://vernier.com/spectral-analysis)

### Go Direct Tris-Compatible Flat pH

This pH sensor features a sealed, gel-filled, double-junction electrode, making it compatible with Tris buffers and solutions containing proteins or sulfides.

GDX-FPH \$139

[vernier.com/gdx-fph](http://vernier.com/gdx-fph)



### BlueView Transilluminator

BlueView™ Transilluminator uses super-bright blue LEDs to illuminate electrophoresis gels stained with fluorescent dyes (e.g., SYBR® Safe). This combination is a safer alternative to ethidium bromide and a UV transilluminator.

BLUE-VIEW ⚡ \$425

[vernier.com/blue-view](http://vernier.com/blue-view)



### Go Direct Temperature

This rugged probe measures the temperature of a variety of substances including air, soil, and water.

Range: -40 to 125°C

GDX-TMP \$79

[vernier.com/gdx-tmp](http://vernier.com/gdx-tmp)



### Go Direct Conductivity

Use this sensor to measure total dissolved solids (TDS) in aquatic samples or the salinity of soil samples.

GDX-CON \$129

[vernier.com/gdx-con](http://vernier.com/gdx-con)



## Vernier and Bio-Rad



Bio-Rad® combines high-quality supplies, equipment, and curricula with outstanding customer service and technical support—things we believe are important to teachers. Vernier and Bio-Rad enhance classroom experiences with joint experiments and curricula for biotechnology.

Download free sample experiments at [vernier.com/bio-rad-kits](http://vernier.com/bio-rad-kits)

# Environmental Science

## Go Direct Optical Dissolved Oxygen

Use this sensor to measure dissolved oxygen, water temperature, and atmospheric pressure. It is ideal for experiments in environmental science.

GDX-ODO \$385

[vernier.com/gdx-odo](http://vernier.com/gdx-odo)



## Go Direct Conductivity

Use this sensor to measure total dissolved solids (TDS) in aquatic samples or the salinity of soil samples.

GDX-CON \$129

[vernier.com/gdx-con](http://vernier.com/gdx-con)



## Go Direct Nitrate Ion-Selective Electrode

Use this sensor to measure nitrate concentration in water samples from water sources throughout your watershed.

GDX-N03 \$299

[vernier.com/gdx-no3](http://vernier.com/gdx-no3)



## Go Direct Tris-Compatible Flat pH

The flat glass shape of this pH sensor is more durable and easier to clean than the traditional pH bulb shape, making it the best choice for environmental science.

GDX-FPH \$139

[vernier.com/gdx-fph](http://vernier.com/gdx-fph)



## Go Direct Salinity

Easily and precisely measure the total dissolved salt content in an aqueous solution. Measure water with a wide variety of salinities, from brackish water to ocean water, and even hyper-saline environments.

GDX-SAL \$180

[vernier.com/gdx-sal](http://vernier.com/gdx-sal)



## Go Direct Soil Moisture

This sensor uses capacitance to measure the volumetric water content of soil.

GDX-SM \$149

[vernier.com/gdx-sm](http://vernier.com/gdx-sm)



## Go Direct Weather System

Easily monitor a wide variety of environmental factors with just one sensor. Go Direct® Weather System includes an affordable, wireless handheld sensor that measures ambient temperature, humidity, wind speed, wind chill, dew point, barometric pressure, and more.

GDX-WTVA \$164

[vernier.com/gdx-wtva](http://vernier.com/gdx-wtva)



## Go Direct PAR

Go Direct PAR (Photosynthetically Active Radiation) measures photosynthetic light levels in both air and water.

GDX-PAR \$299

[vernier.com/gdx-par](http://vernier.com/gdx-par)



## Environmental Science LabQuest 3 Package

This package includes LabQuest® 3 with 10 sensors that work with Vernier Graphical Analysis® Pro and LabQuest 3.

- LabQuest 3
- Stainless Steel Temperature Probe
- Go Direct Tris-Compatible Flat pH
- Go Direct Conductivity Probe
- Go Direct Optical Dissolved Oxygen Probe
- Go Direct Soil Moisture Sensor
- Turbidity Sensor
- Go Direct Weather System
- Go Direct Light and Color Sensor
- Go Direct CO<sub>2</sub> Gas Sensor
- Go Direct Energy Sensor

LQ3-EV-DX \$1,932

[vernier.com/lq3-ev-dx](http://vernier.com/lq3-ev-dx)

Buy 8 or more packages at \$1,894 and save \$304.



## Go Direct Sensor Clamp

Prevent accidental drops during field investigations with the Go Direct Sensor Clamp.

GDX-CLAMP \$14


[vernier.com/gdx-clamp](http://vernier.com/gdx-clamp)



## Renewable Energy with Vernier

The *Renewable Energy with Vernier* lab book features 26 experiments in wind and solar energy. The book contains a combination of explorations, classic experiments, inquiry investigations, engineering projects, and more.

[vernier.com/rev](http://vernier.com/rev)

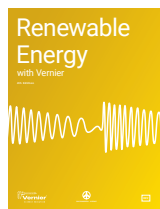

**Printed book + download**  
REV \$59

**Download only**  
REV-E \$45

## Investigating Environmental Science through Inquiry

This book contains 34 inquiry-based environmental science investigations.\* Topics include Earth systems and resources, the living world, global change and population, energy resources and consumption, and pollution.

[vernier.com/esi](http://vernier.com/esi)




**Printed book + download**  
ESI \$59

**Download only**  
ESI-E \$45

## Water Quality with Vernier

With the 18 water quality tests in *Water Quality with Vernier*,\* students investigate the water quality of a body of water by testing pH, total dissolved solids, dissolved oxygen, BOD, and more.

[vernier.com/wqw](http://vernier.com/wqw)

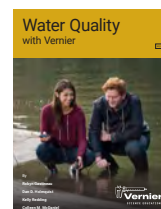

**Printed book + download**  
WQV \$59

**Download only**  
WQV-E \$45

\* Instructions for Vernier Graphical Analysis app and Go Direct sensors are not yet available.

# Biology Products

## Go Direct Sensors

Product	Order Code	Price
Go Direct® Blood Pressure	GDX-BP	\$125
Go Direct CO <sub>2</sub> Gas	GDX-CO2	\$229
Go Direct Colorimeter	GDX-COL	\$135
Go Direct Conductivity	GDX-CON	\$129
Go Direct EKG	GDX-EKG	\$175
Go Direct Ethanol Vapor	GDX-ETOH	\$169
Go Direct Energy	GDX-NRG	\$99
Go Direct Force and Acceleration	GDX-FOR	\$119
Go Direct Gas Pressure	GDX-GP	\$105
Go Direct Hand Dynamometer	GDX-HD	\$135
Heart Rate Monitors		
Go Wireless Exercise Heart Rate	GW-EHR	\$79
Go Wireless® Heart Rate	GW-HR	\$99
Ion-Selective Electrodes		
Go Direct Ammonium Ion-Selective Electrode	GDX-NH4	\$299
Go Direct Nitrate Ion-Selective Electrode	GDX-NO3	\$299
Go Direct Light and Color	GDX-LC	\$99
Go Direct O <sub>2</sub> Gas	GDX-O2	\$219
Go Direct Optical Dissolved Oxygen	GDX-ODO	\$385
Go Direct PAR	GDX-PAR	\$299
pH Sensors		
Go Direct pH	GDX-PH	\$109
Go Direct Tris-Compatible Flat pH	GDX-FPH	\$139
Go Direct Pyranometer	GDX-PYR	\$299
Go Direct Respiration Belt	GDX-RB	\$125
Go Direct Salinity	GDX-SAL	\$180
Go Direct Soil Moisture	GDX-SM	\$149
Spectrophotometers		
Go Direct SpectroVis® Plus	GDX-SVISPL	\$479
Go Direct Fluorescence/ UV-VIS Spectrophotometer	GDX-SPEC-FUV	\$3,499
Go Direct UV-VIS Spectrophotometer	GDX-SPEC-UV	\$2,999
Go Direct Spirometer	GDX-SPR	\$229
Temperature Probes		
Go Direct Surface Temperature	GDX-ST	\$99
Go Direct Temperature	GDX-TMP	\$79
Go Direct Wide-Range Temperature	GDX-WRT	\$139
Go Direct Weather	GDX-WTHR	\$135
Go Direct Weather System	GDX-WTVA	\$164

## LabQuest Sensors

Product	Order Code	Price
PAR Sensor	PAR-BTA	\$259
Relative Humidity Sensor	RH-BTA	\$79
Salinity Sensor	SAL-BTA	\$149
Soil Moisture Sensor	SMS-BTA	\$148
Turbidity Sensor	TRB-BTA	\$124

## Accessories and Lab Equipment

Product	Order Code	Price
BioChamber 250	BC-250 †	\$9
BioChamber 2000	BC-2000 †	\$22
BlueView Transilluminator	BLUE-VIEW †	\$425
Disposable Bacteria Filters (pkg. of 10)	SPR-FIL10	\$45
Disposable Mouthpieces (pkg. of 30)	SPR-MP30	\$17
EKG Electrodes (pkg. of 100)	ELEC	\$19
Go Direct Charge Station	GDX-CRG	\$89
Go Direct Sensor Clamp	GDX-CLAMP	\$14
Nose Clip (pkg. of 10)	SPR-NOSE10	\$10
OHAUS® Balances †	vernier.com/ohaus	
Primary Productivity Kit	PPK	\$49
Reflex Hammer Accessory Kit	RFX-ACC	\$29
Stir Station	STIR	\$149
Water Depth Sampler	WDS	\$114
Water Quality Bottles	WQ-BOT	\$28

## Lab Books\*

Product	Order Code	Price
<i>Biology with Vernier</i>	BWV	\$59
<i>Investigating Biology through Inquiry</i>	BIO-I	\$59
<i>Advanced Biology with Vernier</i> (LabQuest® sensors only)	BIO-A	\$59
<i>Human Physiology Experiments: Volume 1</i>	HSB-HP	\$45
<i>Human Physiology Experiments: Volume 2</i>	ALB-HP2	\$45
<i>Investigating Environmental Science through Inquiry</i> (LabQuest sensors only)	ESI	\$59
<i>Renewable Energy with Vernier</i>	REV	\$59
<i>Water Quality with Vernier</i> (LabQuest sensors only)	WQV	\$59

\* Includes printed book and download; also available as a download only

See all our products for biology at [vernier.com/biology](https://www.vernier.com/biology)



# Chemistry

## Why Vernier?

When you teach with Vernier, you're teaching with a complete chemistry solution. From titrations to spectroscopy, our sensors and instrumentation are backed by powerful analytical software, college-level experiments, and unparalleled support.

### Quality

Durable hardware for lab and field use

### Affordable

Designed for education and education budgets

### Versatile

Supports a variety of devices and experiments



The use of these technologies helps to build students' proficiency using instrumentation while providing them with hands-on experience that will better prepare them for careers in the chemistry field.

*Seth Barrett, PhD  
Muskingham University*

[vernier.com/chemistry](https://www.vernier.com/chemistry)



# A Guide to Vernier Data Collection

We're here to support you as an educator as you incorporate data-collection technology into your instruction. See how our products provide you with affordable laboratory solutions designed for student success.

Our Guarantee: Most of our products are protected by a 5-year limited warranty. And after 5 years? We'll make every attempt to repair your equipment.

## What You Need to Get Started with Go Direct Sensors

### Go Direct Sensor

These versatile sensors connect to your device via Bluetooth® wireless technology or USB.

### Device

Go Direct® sensors connect to a wide variety of commonly used devices, including Chromebooks, computers, tablets, smartphones, and LabQuest® 3.

### Software

Vernier Graphical Analysis® Pro  
Vernier Spectral Analysis®  
Vernier Instrumental Analysis®

### Lab Book

Our popular, award-winning lab books provide hundreds of well-tested, customizable experiments. Our lab books come with a generous site license—purchase once and share files across your department.

## What You Need to Get Started with LabQuest 3

### Sensor

#### Go Direct Sensor

These versatile sensors connect to LabQuest 3 via Bluetooth wireless technology or USB.

*Go Direct spectrometers connect only via USB.*

#### LabQuest Sensor

LabQuest sensors connect directly to LabQuest 3 sensor ports (BTA/BTD).

### LabQuest 3

LabQuest 3 serves as a standalone data-collection platform that works with all Vernier sensors.

### Software

LabQuest App (Included)

### Lab Book

Our popular, award-winning lab books provide hundreds of well-tested, customizable experiments. Our lab books come with a generous site license—purchase once and share files across your department.

## Software

### Vernier Graphical Analysis Pro


Our award-winning app enables advanced data analysis, provides opportunities to reinforce and extend learning, and includes features that support hybrid teaching models. We're always working to add and improve features—check out the "What's New" tab in the app.

**1-Year Site License (unlimited seats)** GAP-1YR \$149

**3-Year Site License (unlimited seats)** GAP-3YR \$299

#### Free Trial for Educators

Try out Graphical Analysis Pro for free for 30 days. Access the sample experiments and enhanced analysis tools to use with your students. Get a free trial and learn about site license options at [vernier.com/graphical-analysis](http://vernier.com/graphical-analysis)

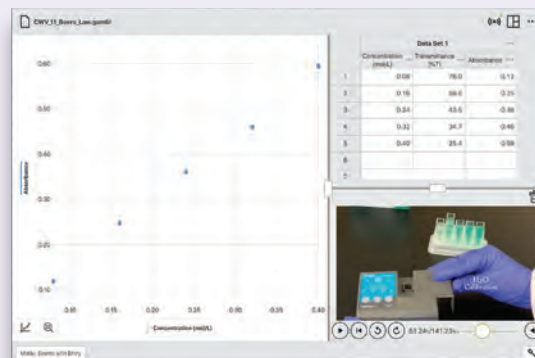

### Vernier Instrumental Analysis

Vernier Instrumental Analysis is used for more advanced instrumentation such as Go Direct Mini GC,™ Go Direct Polarimeter, and Go Direct Cyclic Voltammetry System.

[vernier.com/instrumental-analysis](http://vernier.com/instrumental-analysis)


### Vernier Spectral Analysis

Vernier Spectral Analysis supports our family of spectrometers. Use it to generate full spectra, conduct Beer's law investigations, and investigate kinetics.

[vernier.com/spectral-analysis](http://vernier.com/spectral-analysis)

Get everything you need with our new 4-in-1 software bundle!  
**3-Year Site License** DCA-3YR \$499 | Learn more at [vernier.com/dca-3yr](http://vernier.com/dca-3yr)

## LabQuest 3

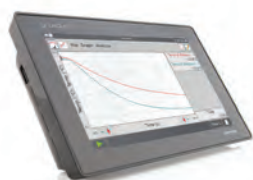

**LabQuest 3 is a powerful, advanced, easy-to-navigate, and versatile data-logging solution for STEM students.**

LabQuest 3 is a standalone data-collection platform that students can use to collect, analyze, and interact with data efficiently. With its touch-screen abilities, students can navigate the platform with ease, and because of its wireless capabilities, students can collect data anywhere.

LABQ3 \$399 [vernier.com/labq3](http://vernier.com/labq3)

### LabQuest App

LabQuest 3 has built-in software that gives your students real-time graphing capabilities in a handheld device. It's powerful, yet beautifully simple.

## General Chemistry

[vernier.com/general-chemistry](http://vernier.com/general-chemistry)

### Go Direct Temperature

Use this rugged temperature probe for investigating endothermic and exothermic reactions, determining the physical properties of water, and investigating intermolecular forces.

Range: -40 to 125°C

GDX-TMP \$79

[vernier.com/gdx-tmp](http://vernier.com/gdx-tmp)



### Go Direct pH

Go Direct pH is an important and versatile sensor for your laboratory. Conduct acid-base titrations, monitor pH changes during chemical reactions, and investigate buffers. The wireless connection makes it easier to do field-based studies such as testing the pH of surface water.

GDX-PH \$109

[vernier.com/gdx-ph](http://vernier.com/gdx-ph)



### Go Direct Gas Pressure

Explore gas laws and the Clausius-Clapeyron equation with this sensor that measures the absolute pressure of a gas.

Range: 0 to 400 kPa

GDX-GP \$105

[vernier.com/gdx-gp](http://vernier.com/gdx-gp)



### Go Direct SpectroVis® Plus

With a range of 380 to 950 nm, students can use this spectrophotometer to easily collect a full-wavelength spectrum, study absorbance vs. concentration, or monitor rates of reaction.

Collect and analyze data using Vernier Spectral Analysis or LabQuest App.

GDX-SVISPL \$479

[vernier.com/gdx-svispl](http://vernier.com/gdx-svispl)



### Go Direct Drop Counter

As an alternative to using a buret, the drop counter precisely records the number of drops of titrant added during a titration and then automatically converts it to volume.

GDX-DC \$115

[vernier.com/gdx-dc](http://vernier.com/gdx-dc)



### Stir Station

This combination stir plate/ring stand can be used with AC power (included) or four C batteries (not included).

STIR \$149

[vernier.com/stir](http://vernier.com/stir)




### Advanced Chemistry with Vernier

This book contains 35 ready-to-use student experiments that support general chemistry. Instructor notes with sample data are also included.

Topics

- Gas laws
- Titrations
- Spectroscopy
- Electrochemistry

[vernier.com/chem-a](http://vernier.com/chem-a)

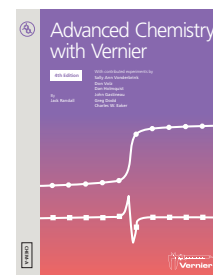

Printed book + download

CHEM-A \$59

Download only


CHEM-A-E \$45

## Go Direct Polarimeter

The concept of chirality can be difficult for students to visualize. Go Direct® Polarimeter provides a visual representation of this concept by measuring the optical rotation of optical isomers such as sugars, amino acids, and proteins.

GDX-POL \$545

Free experiment downloads available at [vernier.com/gdx-pol](http://vernier.com/gdx-pol)

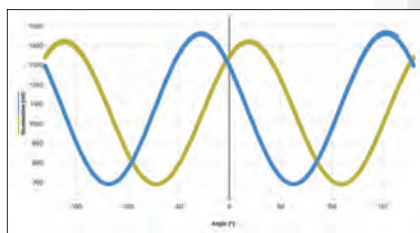

Comparing the optical rotations of fructose and sucrose

### Free Software

Vernier Instrumental Analysis

See page 15.

## Go Direct Fluorescence/UV-VIS Spectrophotometer

This spectrophotometer measures the fluorescence and absorbance spectra of samples such as quinine sulfate, fluorescein, chlorophyll, and fluorescent proteins, all while connecting to your device via Bluetooth® wireless technology or USB.

GDX-SPEC-FUV \$3,499

Free experiment downloads available at

[vernier.com/gdx-spec-fuv](http://vernier.com/gdx-spec-fuv)

### Wavelength Range

- 220 to 850 nm

### Light Sources

- Visible: LED-boosted tungsten
- UV: Deuterium
- Fluorescence: Exchangeable LEDs for excitation at 375 nm, 450 nm, and 525 nm (additional wavelengths sold separately)




## Go Direct Tris-Compatible Flat pH

Go Direct Tris-Compatible Flat pH is a double-junction electrode for measuring pH in Tris buffers and solutions containing proteins or sulfides. The flat glass shape makes it easy to clean and useful for measuring the pH of semisolids such as soil slurries and certain foods.

GDX-FPH \$139

[vernier.com/gdx-fph](http://vernier.com/gdx-fph)

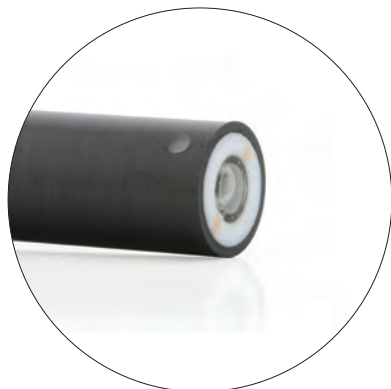


## Vernier Spectral Analysis


Our free Vernier Spectral Analysis® app makes it easy to incorporate spectroscopy into your chemistry experiments. Using the app, students can collect a full spectrum and explore topics such as Beer's law, kinetics, and fluorescence.

The user-friendly software includes analysis features such as curve fitting and data interpolation.

[vernier.com/spectral-analysis](http://vernier.com/spectral-analysis)

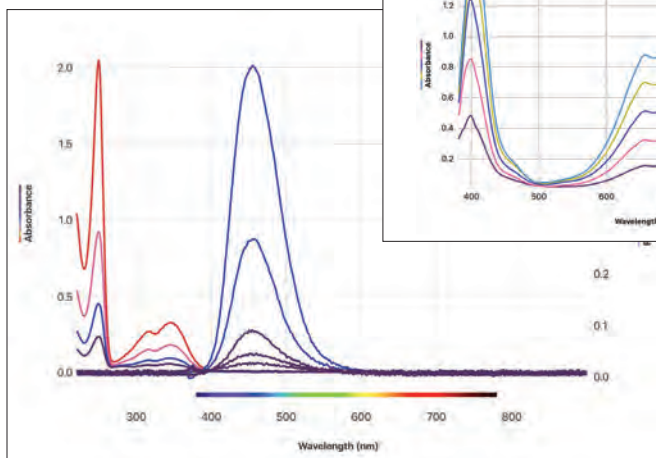

Quinine sulfate spectra at varying concentrations; absorbance (left) and fluorescence with excitation at 375 nm (right)

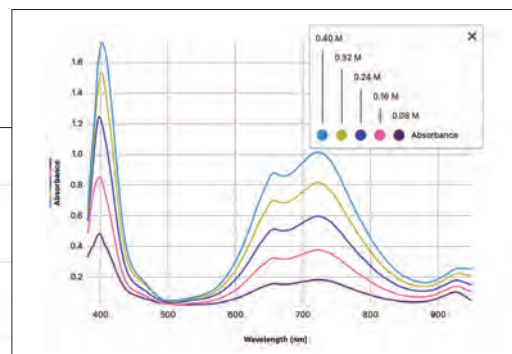

Nickel sulfate absorbance spectra at varying concentrations

## Go Direct Melt Station

Go Direct Melt Station accurately measures melting temperatures of a solid (up to 260°C), and real-time graphing provides a unique perspective of the melting process.

GDX-MLT  \$599

[vernier.com/gdx-mlt](http://vernier.com/gdx-mlt)



## Go Direct UV-VIS Spectrophotometer

The Go Direct UV-VIS Spectrophotometer connects to your device via Bluetooth wireless technology or USB to generate full spectra, Beer's law data, and kinetic traces of ultraviolet and visible-absorbing samples such as aspirin, DNA, proteins, and NADH.

GDX-SPEC-UV \$2,999

[vernier.com/gdx-spec-uv](http://vernier.com/gdx-spec-uv)

### Wavelength Range

- 220 to 850 nm

### Light Sources

- Visible: LED-boosted tungsten
- UV: Deuterium



## Go Direct Wide-Range Temperature

Go Direct Wide-Range Temperature is designed to be used as you would use a thermometer for experiments such as the recrystallization of benzoic acid, simple and fractional distillations, determination of boiling points, the synthesis and analysis of aspirin and other organic compounds, and more.

Range: -20 to 330°C

GDX-WRT \$139

[vernier.com/gdx-wrt](http://vernier.com/gdx-wrt)



### Free Software

#### Vernier Spectral Analysis

See page 14.

## Go Direct Mini GC

With the easy-to-use Go Direct Mini GC™ and the free Vernier Instrumental Analysis® app, students can separate, analyze, and identify substances contained in a volatile liquid or gaseous sample. This portable gas chromatograph detects polar and nonpolar compounds allowing for a wide range of experiments. Sample experiments include fractional distillation and Fischer esterification.

GDX-GC \$3,999

[vernier.com/gdx-gc](http://vernier.com/gdx-gc)




## Organic Chemistry with Vernier

*Organic Chemistry with Vernier* contains 26 experiments that represent a broad range of topics and techniques taught in most college organic chemistry lab courses. The experiments in this book build upon prior knowledge, laboratory techniques, and skills that students learn in general chemistry courses.

### Topics

- Distillation
- Chromatography
- Synthesis
- Polarimetry

[vernier.com/chem-o](http://vernier.com/chem-o)

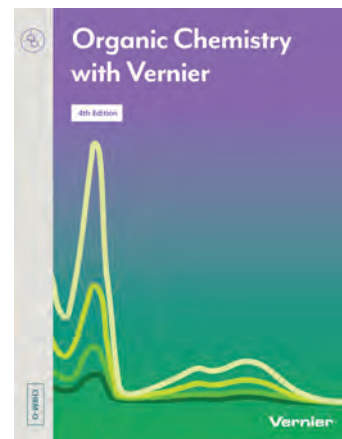

### Printed book + download


CHEM-0 \$59

### Download only

CHEM-0-E \$45

## Vernier Instrumental Analysis

With our free Vernier Instrumental Analysis app, students can collect and analyze data from our Go Direct Mini GC, Go Direct Polarimeter, and Go Direct Cyclic Voltammetry System (page 16) using computers, Chromebooks, or compatible mobile devices.

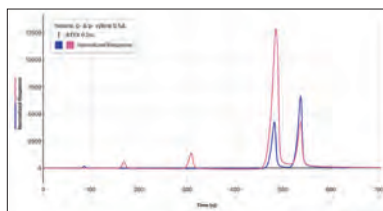

Comparing chromatograms of a BTEX mixture with hexane and xylene isomers

[vernier.com/instrumental-analysis](http://vernier.com/instrumental-analysis)



# Analytical and Inorganic Chemistry

## Go Direct Cyclic Voltammetry System

Give your students hands-on experience with electrochemically active reactions using this affordable potentiostat and disposable screen-printed electrodes.

Easily incorporate electrochemistry into your curriculum using our e-book, *Electrochemistry Experiments with the Go Direct® Cyclic Voltammetry System*, available for free with your purchase.

GDX-CVS \$999

Free experiment downloads available at [vernier.com/gdx-cvs](http://vernier.com/gdx-cvs)



### Free Software

**Vernier Instrumental Analysis®**

See page 15.

## NEW Accessories for an Eco-Friendly Lab

Three new products are now a part of the Vernier college chemistry lineup, offering reliable, easy-to-use solutions that enhance the lab experience for students and educators. These accessories can be used alongside Vernier Go Direct sensors to provide impactful learning experiences that are both efficient and environmentally conscious.

### Electrochemistry Half-Cell Plate

The Electrochemistry Half-Cell Plate is a reusable, low-volume platform that holds up to eight metal-ion half-cells, each accessing the central salt bridge. When paired with the Electrochemistry Metals Kit and Go Direct Voltage Probe, it ensures precise, consistent exploration of oxidation-reduction reactions, the electrochemical series, and Nernst variables, while simplifying the setup and cleanup of electrochemical experiments.

ECHEM-PLT \$59

[vernier.com/echem-plt](http://vernier.com/echem-plt)

### Electrochemistry Metals Kit

The Electrochemistry Metals Kit provides a convenient, reusable set of seven different metal samples for conducting electrochemical experiments (aluminum, copper, iron, lead, nickel, silver, and zinc). When paired with a half-cell plate and voltage probe, this convenient kit simplifies the study of oxidation, reduction, the electrochemical series, and Nernst variables.

ECHEM-MTLS ⬆️ \$20



[vernier.com/echem-mlts](http://vernier.com/echem-mlts)

### Precision Volume Dispenser

The Precision Volume Dispenser gives students and educators a precise, rugged, and easy-to-use alternative to a burette. It can be used with the Go Direct Drop Counter, Go Direct pH Sensor, Go Direct Conductivity Probe, and Go Direct Oxidation-Reduction Potential (ORP) Sensor to conduct acid-base, conductometric, or potentiometric titrations.

DC-DISP \$85

[vernier.com/dc-disp](http://vernier.com/dc-disp)

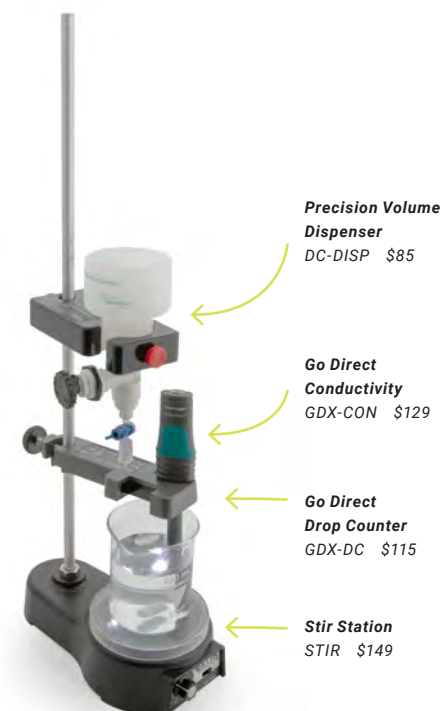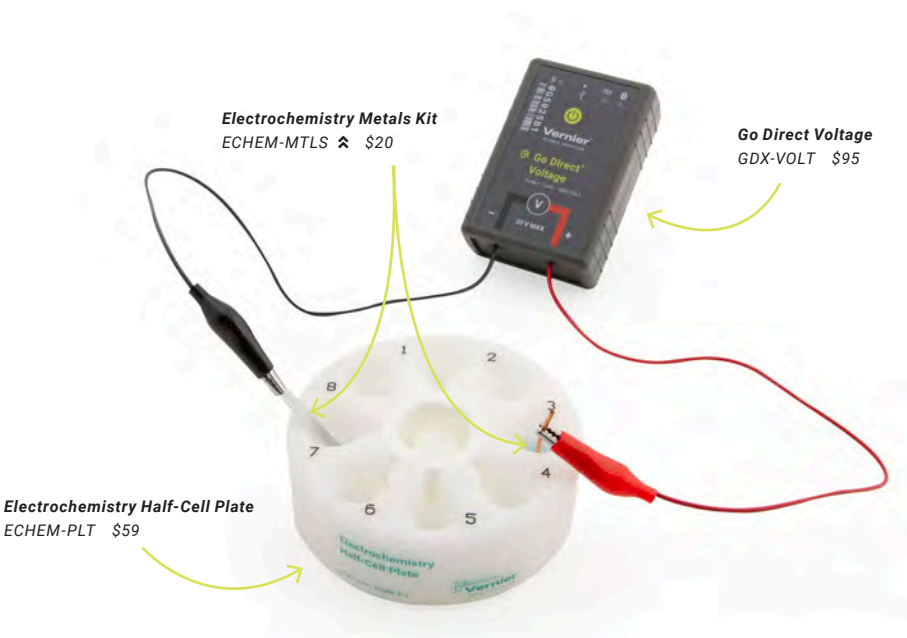



## Go Direct Fluorescence/UV-VIS Spectrophotometer

The Go Direct Fluorescence/UV-VIS Spectrophotometer measures the fluorescence and absorbance spectra of various samples.

GDX-SPEC-FUV \$3,499

### Wavelength Range

- 220 to 850 nm

### Light Sources

- Visible: LED-boosted tungsten
- UV: Deuterium
- Fluorescence: Exchangeable LEDs for excitation at 375 nm, 450 nm, and 525 nm (additional wavelengths sold separately)

Free experiment downloads available at [vernier.com/gdx-spec-fuv](http://vernier.com/gdx-spec-fuv)



### Free Software

Vernier Spectral Analysis®

See page 14.

## Vernier Flash Photolysis Spectrometer

The Vernier Flash Photolysis Spectrometer is perfect for students to explore the fundamental principles of photochemical reactions.

VSP-FP \$5,999

### Wavelength Range

- 450 to 750 nm

### Light Sources

- Xenon flashlamp (pump) white LED (probe)
- 13 exchangeable filters for detection wavelength
- 2 exchangeable colored glass filters for excitation wavelength

Free experiment downloads available at [vernier.com/vsp-fp](http://vernier.com/vsp-fp)



### Free Software

Collect data with the free, all-inclusive Flash Photolysis Spectrometer Software.\*

\*Available for Windows® only

## Go Direct Visible Spectrophotometer

Conduct full-spectrum, equilibrium, kinetic, concentration absorbance studies, or emissions studies with the Go Direct Visible Spectrophotometer. Its quality optical bench delivers high-accuracy measurements, and the aluminum housing minimizes temperature fluctuations.

GDX-SPEC-VIS \$1,999

### Wavelength Range

- 380 to 950 nm

### Light Sources

- Visible: LED-boosted tungsten

[vernier.com/gdx-spec-vis](http://vernier.com/gdx-spec-vis)



### Free Software

Vernier Spectral Analysis

See page 14.

## Go Direct Mini GC

With the easy-to-use Go Direct Mini GC™ and the free Vernier Instrumental Analysis app, students can separate, analyze, and identify substances contained in a volatile liquid sample.

GDX-GC \$3,999

[vernier.com/gdx-gc](http://vernier.com/gdx-gc)



## Go Direct Wide-Range Pressure

Measure absolute pressures up to 690 kPa (100 psi) with excellent accuracy. The Gas and Liquid Compatible Go Direct Wide-Range Pressure (order code GDX-WRPL) can measure changes in liquid, gas, or vapor pressure. This liquid-compatible version is an excellent choice for agricultural and environmental studies, particularly those involving hydraulic lift and fluid pressure measurements. Go Direct Wide-Range Pressure (order code GDX-WRP) can measure changes in gas or vapor pressure.



Liquid & Gas Compatible GDX-WRPL \$249  
[vernier.com/gdx-wrpl](http://vernier.com/gdx-wrpl)



Gas Compatible GDX-WRP \$219  
[vernier.com/gdx-wrp](http://vernier.com/gdx-wrp)

# Chemistry Products

## Go Direct Sensors

Sensor	Order Code	Price
Go Direct® Colorimeter	GDX-COL	\$135
Conductivity Probes		
Go Direct Conductivity	GDX-CON	\$129
Go Direct Platinum-Cell Conductivity	GDX-CONPT	\$189
Current Probes		
Go Direct Constant Current System	GDX-CCS	\$95
Go Direct Current	GDX-CUR	\$99
Go Direct Drop Counter	GDX-DC	\$115
Gas Pressure Sensors		
Go Direct Gas Pressure	GDX-GP	\$105
Go Direct Wide-Range Pressure (Gas Compatible)	GDX-WRP	\$219
Go Direct Wide-Range Pressure (Liquid & Gas Compatible)	GDX-WRPL	\$249
Go Direct ORP	GDX-ORP	\$119
pH Sensors		
Go Direct Glass-Body pH	GDX-GPH	\$159
Go Direct pH	GDX-PH	\$109
Go Direct Tris-Compatible Flat pH	GDX-FPH	\$139
Go Direct Radiation Monitor	GDX-RAD	\$219
Temperature Probes		
Go Direct Surface Temperature	GDX-ST	\$99
Go Direct Temperature	GDX-TMP	\$79
Go Direct Thermocouple	GDX-TC	\$118
Go Direct Wide-Range Temperature	GDX-WRT	\$139
Go Direct Voltage	GDX-VOLT	\$95

## LabQuest Sensors

Sensor	Order Code	Price
Colorimeter	COL-BTA	\$128
Conductivity Probe	CON-BTA	\$129
Current Probes		
Constant Current System	CCS-BTA	\$79
Current Probe	DCP-BTA	\$49
Drop Counter	VDC-BTD	\$109
Electrode Amplifier	EA-BTA	\$55
Gas Pressure Sensors		
Gas Pressure Sensor	GPS-BTA	\$103
Pressure Sensor 400	PS400-BTA	\$229
Instrumentation Amplifier	INA-BTA	\$89
Melt Station	MLT-BTA ⚡	\$629
ORP Sensor	ORP-BTA	\$115
pH Sensors		
Glass-Body pH Electrode BNC (requires Electrode Amplifier)	GPH-BNC	\$95
pH Sensor	PH-BTA	\$105
Tris-Compatible Flat pH Sensor	FPH-BTA	\$124
Radiation Monitor	VRM-BTD	\$199
Temperature Probes		
Stainless Steel Temperature Probe	TMP-BTA	\$45
Surface Temperature Sensor	STS-BTA	\$34
Thermocouple	TCA-BTA	\$79
Wide-Range Temperature Probe	WRT-BTA	\$119
Voltage Probes		
Differential Voltage Probe	DVP-BTA	\$49
Voltage Probe	VP-BTA	\$14

## Instrumentation

Instrument	Order Code	Price
Go Direct Cyclic Voltammetry System	GDX-CVS	\$999
Go Direct Melt Station	GDX-MLT ⚡	\$599
Go Direct Mini GC™	GDX-GC	\$3,999
Go Direct Polarimeter	GDX-POL	\$545

## Spectrometers

[vernier.com/spectrometers](http://vernier.com/spectrometers)

Spectrometer	Order Code	Price
Go Direct Emissions Spectrometer	GDX-SPEC-EM	\$950
Vernier Flash Photolysis Spectrometer	VSP-FP	\$5,999
Go Direct Fluorescence/UV-VIS Spectrophotometer	GDX-SPEC-FUV	\$3,499
Go Direct SpectroVis® Plus	GDX-SVISPL	\$479
Go Direct UV-VIS Spectrophotometer	GDX-SPEC-UV	\$2,999
Go Direct Visible Spectrophotometer	GDX-SPEC-VIS	\$1,999

## Lab Equipment

Equipment	Order Code	Price
Electrode Support	ESUP	\$10
OHAUS® Balances ⚡	<a href="http://vernier.com/ohaus">vernier.com/ohaus</a>	
Stir Station	STIR	\$149
Cuvette Rack	CUV-RACK	\$9
Go Direct Charge Station	GDX-CRG	\$89
<b>NEW</b> Precision Volume Dispenser	DC-DISP	\$85
<b>NEW</b> Electrochemistry Half-Cell Plate	ECHEM-PLT	\$59
<b>NEW</b> Electrochemistry Metals Kit	ECHEM-MTLS ⚡	\$20

## Lab Books\*

Product	Order Code	Price
<i>Advanced Chemistry with Vernier</i>	CHEM-A	\$59
<i>Chemistry with Vernier</i>	CWV	\$59
<i>Food Chemistry Experiments</i>	HSB-FOOD	\$45
<i>Forensic Chemistry Experiments</i>	HSB-FCHEM	\$45
<i>Investigating Chemistry through Inquiry</i>	CHEM-I	\$59
<i>Organic Chemistry with Vernier</i>	CHEM-O	\$59

\* Includes printed book and download; also available as a download only

See all our products for college chemistry online at [vernier.com/college-chemistry](http://vernier.com/college-chemistry)

# Physics

## Why Vernier?

Vernier started when one educator, Dave Vernier, decided to build solutions to bring physics to life for his students. Today, our complete physics solution is still powered by the desire to inspire students and foster learning and is backed by powerful software and unparalleled support.

### Quality

Durable hardware lasts for years of use

### Affordable

Designed for education and education budgets

### Versatile

Supports a variety of devices and experiments



“

I find your extensive suite of hardware—and your innovative software—extremely helpful in my teaching. Couldn't do it without your stuff!

*Barbara Hughey  
Massachusetts Institute of Technology  
Cambridge, MA*

[vernier.com/physics](http://vernier.com/physics)

# A Guide to Vernier Data Collection

We're here to support you as an educator as you incorporate data-collection technology into your instruction. See how our products provide you with affordable laboratory solutions designed for student success.

Our Guarantee: Most of our products are protected by a 5-year limited warranty. And after 5 years? We'll make every attempt to repair your equipment.

## What You Need to Get Started with Go Direct Sensors

### Go Direct Sensors

These versatile sensors connect to your device via Bluetooth® wireless technology or USB.

### Device

Go Direct® sensors connect to a wide variety of commonly used devices, including Chromebooks, computers, smartphones, tablets, and LabQuest® 3.

### Software

Vernier Graphical Analysis® Pro  
Vernier Spectral Analysis®  
Vernier Video Analysis®

### Lab Book

Our popular, award-winning lab books provide hundreds of well-tested, customizable experiments. Our lab books come with a generous site license—purchase once and share files across your department.

## What You Need to Get Started with LabQuest Sensors

### LabQuest Sensor

LabQuest sensors share data with your device via a wired connection (BTA/BTD) to an interface from the LabQuest family.

### Interface

An interface sends information from the sensor to the data-collection and analysis software. The LabQuest family includes LabQuest 3 and LabQuest Mini.

### Software

Vernier Graphical Analysis Pro  
LabQuest App  
Vernier Video Analysis

### Lab Book

Our popular, award-winning lab books provide hundreds of well-tested, customizable experiments.

### Device

LabQuest sensors connect to computers, Chromebooks, and compatible mobile devices through a LabQuest interface.

## Software

### Vernier Graphical Analysis Pro


Our award-winning app enables advanced data analysis, provides opportunities to reinforce and extend learning, and includes features that support hybrid teaching models. We're always working to add and improve features—check out the “What’s New” tab in the app.

**1-Year Site License (unlimited seats)** GAP-1YR \$149

**3-Year Site License (unlimited seats)** GAP-3YR \$299

#### Free Trial for Educators

Try out Graphical Analysis Pro for free for 30 days. Access the sample experiments and enhanced analysis tools to use with your students. Get a free trial and learn about site license options at [vernier.com/graphical-analysis](http://vernier.com/graphical-analysis)

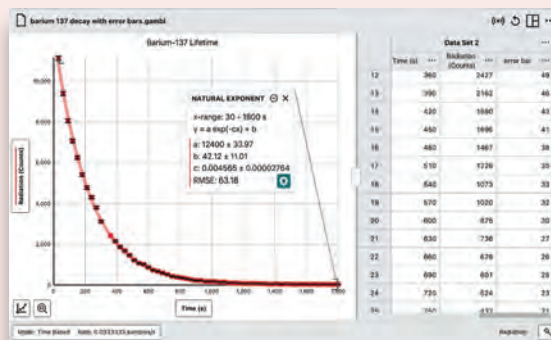

### Vernier Spectral Analysis

Spectral Analysis supports our family of spectrometers. The user-friendly interface walks students through the data-collection process and includes analysis features such as curve fitting and data interpolation.

[vernier.com/spectral-analysis](http://vernier.com/spectral-analysis)

### Vernier Video Analysis

Students use their smartphone or tablet in the laboratory or out in the field to record motion. They then import the video into Video Analysis on any device to mark the object in motion, set the scale, and create graphs of the motion.

[vernier.com/video-analysis](http://vernier.com/video-analysis)

Get everything you need with our new 4-in-1 software bundle!

**3-Year Site License** DCA-3YR \$499 | Learn more at [vernier.com/dca-3yr](http://vernier.com/dca-3yr)

## LabQuest 3

LabQuest® 3 is a standalone data-collection platform that students can use to collect, analyze, and interact with data efficiently. It can also double as a sensor interface for computers.

LABQ3 \$399 [vernier.com/labq3](http://vernier.com/labq3)

## LabQuest Mini

LabQuest Mini is a sensor interface that brings you the power of our award-winning LabQuest technology when you don't need the versatility of a standalone device.

LQ-MINI \$189 [vernier.com/lq-mini](http://vernier.com/lq-mini)





## Free 30-Day Preview

Not sure if our technology is right for your laboratory? Our risk-free preview program makes it easier to decide. Most Vernier products are available for a 30-day preview (or longer, if requested) to educational institutions.

# Mechanics

## Dynamics Cart and Track System with Go Direct Sensor Cart

The Dynamics Cart and Track System (DTS) with Go Direct Sensor Cart includes essential laboratory equipment for teaching dynamics and kinematics. With our Go Direct Sensor Cart, students can explore force, position, velocity, and acceleration directly on their device using Bluetooth wireless technology. There are no wires to create drag, and no additional equipment is required! Each cart features built-in sensors that simplify experiment setup and make this system the best choice for studying dynamics and kinematics.

with 1.2 m Track DTS-GDX \$636  
with 2.2 m Track DTS-GDX-LONG \$748

[vernier.com/dts-gdx](http://vernier.com/dts-gdx)

Additional Cart and Track options are available at [vernier.com/dynamics](http://vernier.com/dynamics)



## Cart Fan

The Cart Fan is a modular fan you can add to DTS and Go Direct carts that provides a constant force for dynamics investigations. Using one or multiple fans on a single cart, students can investigate constant acceleration, balanced/unbalanced forces, and variable thrust angles.

DTS-CFAN \$89

[vernier.com/dts-cfan](http://vernier.com/dts-cfan)



## Go Direct Force and Acceleration

Measure forces as small as  $\pm 0.1$  N and up to  $\pm 50$  N with this sensor that couples a 3-axis accelerometer with a stable and accurate force sensor. It also includes a 3-axis gyroscope for experiments involving rotation.

GDX-FOR \$119

[vernier.com/gdx-for](http://vernier.com/gdx-for)



## Go Direct Centripetal Force System

Conduct a variety of rotational dynamics investigations with a single experiment setup. This combination easily measures angular velocity, centripetal force, and acceleration without the hassle of wires.

GDX-CFAF \$444

[vernier.com/gdx-cfaf](http://vernier.com/gdx-cfaf)



## Moment of Inertia Accessory Kit

With the Moment of Inertia Accessory Kit, students can explore inertia in a broader context. The kit expands the capabilities of the Vernier centripetal force apparatus when investigating moments of inertia of different geometries.

CFA-MIK \$219

[vernier.com/cfa-mik](http://vernier.com/cfa-mik)



## Go Direct Motion

Use ultrasound to measure the position, velocity, and acceleration of moving objects.

GDX-MD \$129

[vernier.com/gdx-md](http://vernier.com/gdx-md)



## Go Direct Force Plate

Measure the forces developed during stepping, jumping, and other human-scale actions. Observe change in apparent weight as you ride an elevator, or measure reaction forces as you lean against a wall.

GDX-FP \$399

[vernier.com/gdx-fp](http://vernier.com/gdx-fp)



## Go Direct Photogate

This double-gate sensor includes two photogates built into the arms of the sensor. It accurately measures velocity and acceleration.

GDX-VPG \$105

[vernier.com/gdx-vpg](http://vernier.com/gdx-vpg)





# Waves and Sound

## Go Direct Sound

This is really two sensors in one—measure sound level in decibels, or capture and evaluate sound waveforms.

GDX-SND \$99

[vernier.com/gdx-snd](http://vernier.com/gdx-snd)



## Resonance Apparatus

The Resonance Apparatus is used for the classic experiment of determining the speed of sound using the principle of resonance in a tube that is closed at one end.

RES-APP \$199

[vernier.com/res-app](http://vernier.com/res-app)



Watch the video!



## Power Amplifier

Use the Power Amplifier as a power supply for DC and AC circuit investigations or to drive devices such as speakers.

PAMP \$249

[vernier.com/pamp](http://vernier.com/pamp)



## Power Amplifier Accessory Speaker

Study mechanical waves on strings and springs.

PAAS-PAMP ⚡ \$155

[vernier.com/paas-pamp](http://vernier.com/paas-pamp)



## Frequency Generator

Use the Frequency Generator to create sine, square, sawtooth, and triangle waves at a wide range of frequencies.

FGEN-PAMP \$189

[vernier.com/fgen-pamp](http://vernier.com/fgen-pamp)



# Thermodynamics

## Go Direct Gas Pressure

This sensor measures the absolute pressure of a gas.

GDX-GP \$105

[vernier.com/gdx-gp](http://vernier.com/gdx-gp)



## Go Direct Temperature

Go Direct® Temperature is a durable, stainless steel temperature sensor for use in liquids or air.

Range:  $-40$  to  $125^{\circ}\text{C}$

GDX-TMP \$79

[vernier.com/gdx-tmp](http://vernier.com/gdx-tmp)



## Go Direct Surface Temperature

An exposed temperature sensor makes this an ideal choice for situations where low thermal mass and extremely rapid response are needed. Use this sensor in air and water only.

Range:  $-25$  to  $125^{\circ}\text{C}$

GDX-ST \$99

[vernier.com/gdx-st](http://vernier.com/gdx-st)



# Electricity and Magnetism

## Go Direct Voltage

This sensor combines a wide input voltage range and high precision, making it an excellent choice for investigations of both AC/DC circuits and electromagnetism.

Ranges:  $\pm 20$  V and  $\pm 1$  V

GDX-VOLT \$95

[vernier.com/gdx-volt](http://vernier.com/gdx-volt)



## Go Direct Current

Measure electric currents in circuits with this versatile sensor.

Ranges:  $\pm 1$  A and  $\pm 0.1$  A

GDX-CUR \$99

[vernier.com/gdx-cur](http://vernier.com/gdx-cur)



## Go Direct Static Charge

Unlike a traditional electroscope, Go Direct Static Charge offers a means for quantitative measurement and analysis of positive and negative charges of objects that would not be possible in a traditional lab—such as measuring the charge on a balloon.

Range:  $\pm 100$  nC

GDX-Q \$109

[vernier.com/gdx-q](http://vernier.com/gdx-q)



## Go Direct 3-Axis Magnetic Field

Determine the magnitude and direction of a magnetic field at any point in space with this 3-axis sensor.

Ranges:  $\pm 5$  mT and  $\pm 130$  mT

GDX-3MG \$99

[vernier.com/gdx-3mg](http://vernier.com/gdx-3mg)



## Electrostatics Kit

When using the Electrostatics Kit with Go Direct Static Charge, students can conduct a range of experiments in electrostatics.

ESK-CRG \$129

[vernier.com/esk-crg](http://vernier.com/esk-crg)



## High-Voltage Electrostatics Kit

Use this kit to investigate the distribution of charge on a sphere, transfer of charge on contact between two spheres, and charging by induction.

HVEK-CRG \$299

[vernier.com/hvek-crg](http://vernier.com/hvek-crg)



## Extech® Digital Power Supply

This power supply provides constant current or constant voltage for physics activities that require DC power.

EXPS \$259

[vernier.com/exps](http://vernier.com/exps)



## Vernier Circuit Board 2

Use this convenient platform to study basic series and parallel circuits as well as RLC circuits. Many components for experimentation are provided, and additional components can be added to expand the capability of this useful board.

VCB2 \$145

[vernier.com/vcb2](http://vernier.com/vcb2)



## Electrostatic High-Voltage Genecon

A great addition to the High Voltage Electrostatics Kit, the Electrostatic High-Voltage Genecon generates both positive and negative charges and reliably creates charge differences in high humidity.

HVEK-GEN \$229

[vernier.com/hvek-gen](http://vernier.com/hvek-gen)



# Light and Optics

## Go Direct Light and Color

This sensor combines visible light, UV, and RGB sensors to measure source emission, transmittance, and reflection of light in the visible light to ultraviolet electromagnetic spectrum.

GDX-LC \$99

[vernier.com/gdx-lc](http://vernier.com/gdx-lc)



## Diffraction Apparatus\*†


This set includes a Red Diffraction Laser, Diffraction Slit System, and Combination Linear Position and High Sensitivity Light Sensor. It does not include the required Combination Track/Optics Bench.

DAK \$669 [vernier.com/dak](http://vernier.com/dak)

### Green Diffraction Laser (optional)

Add this laser to your Diffraction Apparatus to study the effect of wavelength on a diffraction pattern.

GDL-DAK \$250 [vernier.com/gdl-dak](http://vernier.com/gdl-dak)

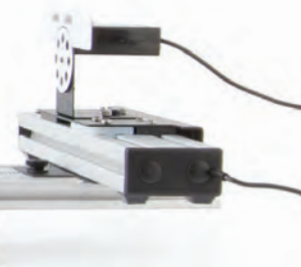

## Optics Expansion Kit

Use the Optics Expansion Kit† with your dynamics track to conduct optics experiments, such as image formation with lenses and light intensity vs. distance. You can even use the kit to build a basic telescope.

OEK \$199

[vernier.com/oek](http://vernier.com/oek)

Kit includes

- 3 lenses (100 mm converging lens, 200 mm converging lens, -150 mm diverging lens)
- Screen
- Power supply
- Combination luminous and point light source
- Light sensor holder
- Aperture screen

The Optics Expansion Kit is used in experiments in our *Physics with Vernier* and *Advanced Physics with Vernier—Beyond Mechanics* lab books.

See website for replacement parts.

\*† Requires an interface such as LabQuest® 3 or LabQuest Mini

† Requires a Combination 1.2 m Track/Optics Bench (TRACK, \$160)

## Accessories

### Color Mixer Kit†

CM-OEK \$189 [vernier.com/cm-oek](http://vernier.com/cm-oek)




### Mirror Set†

M-OEK \$69 [vernier.com/m-oek](http://vernier.com/m-oek)



### Polarizer/Analyzer Set†

PAK-OEK \$89 [vernier.com/pak-oek](http://vernier.com/pak-oek)

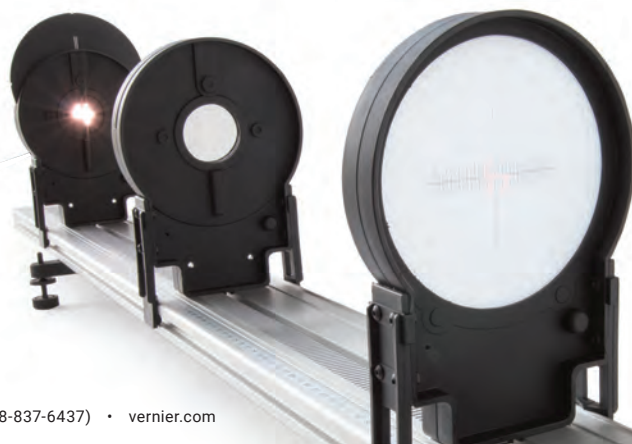



## Go Direct Radiation Monitor

Our radiation monitor detects alpha, beta, gamma, and X-ray radiation. It can be used to explore radiation statistics, measure the rate of nuclear decay, monitor radon progeny, and investigate the effects of shielding. The sensor includes both LED and audible indicators.

GDX-RAD \$219

[vernier.com/gdx-rad](http://vernier.com/gdx-rad)




## Go Direct Emissions Spectrometer

The Go Direct® Emissions Spectrometer connects to your device via Bluetooth® wireless technology or USB to give precise measurements over a range of 350–900 nm. Use it with or without the Vernier Emissions Fiber (not included) to examine spectra of light bulbs, spectrum tubes, or the sun.

GDX-SPEC-EM \$950

[vernier.com/gdx-spec-em](http://vernier.com/gdx-spec-em)

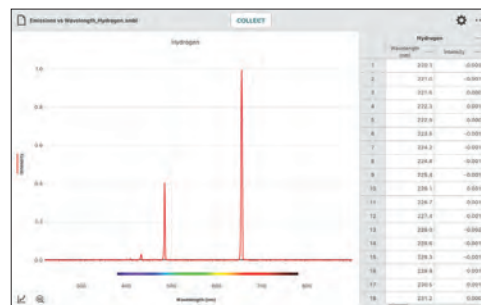

### Vernier Spectral Analysis App

Our free Vernier Spectral Analysis® app combined with our Emissions Spectrometer makes it easy to analyze spectra. Students can quickly locate peaks or compare spectra from different sources.

[vernier.com/spectral-analysis](http://vernier.com/spectral-analysis)



### Vernier Emissions Fiber

VSP-EM-FIBER \$94

[vernier.com/vsp-em-fiber](http://vernier.com/vsp-em-fiber)

## Spectrum Tube Power Supply

This power supply features an ultra-safe design for electrifying spectrum tubes.

ST-SPS ⚡ \$299

[vernier.com/st-sps](http://vernier.com/st-sps)




## Spectrum Tubes

Spectrum Tubes are permanently enclosed in protective plastic carriers, with no exposed high voltage.

Spectrum Tubes carry a warranty of 2 years or 100 hours, whichever comes first (hydrogen tube: 2 years or 40 hours, whichever comes first).

[vernier.com/spectrum-tubes](http://vernier.com/spectrum-tubes)

Spectrum Tube	Order Code	Price
Hydrogen	ST-H ⚡	\$49
Nitrogen	ST-N ⚡	\$49
Helium	ST-HE ⚡	\$49
Neon	ST-NE ⚡	\$49
Carbon Dioxide	ST-CO2 ⚡	\$49
Air	ST-AIR ⚡	\$49
Argon	ST-AR ⚡	\$49

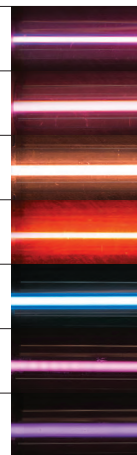

# Additional Physics Products

## Mechanics

Product	Order Code	Price
Go Direct® Acceleration	GDX-ACC	\$115
3-Axis Accelerometer	3D-BTA	\$105
25-g Accelerometer	ACC-BTA	\$110
Bumper and Launcher Kit	BLK	\$99
Cart Fan	DTS-CFAN	\$89
Dual-Range Force Sensor	DFS-BTA	\$124
Dynamics Cart and Track System	DTS	\$360
Dynamics Cart and Track System with Motion Encoder	DTS-EC	\$540
Eddy Current Brake	DTS-ECB	\$19
Encoder Fan Cart	CART-FEC	\$245
Fan Cart	CART-F	\$119
Force Plate	FP-BTA	\$395
Friction Block	PHY-FRBL	\$17
Friction Pad DTS	DTS-PAD	\$32
Go Direct Force Plate	GDX-FP	\$399
Go Direct Sensor Cart Accessory Kit	GDX-CART-AK	\$62
Hanging Mass	PHY-HM250 ⚡	\$23
Independence of Motion	IOM-VPL	\$59
Low-g Accelerometer	LGA-BTA	\$99
Meter Stick	PHY-METER	\$11
Motion Detector	MD-BTD	\$105
Photogate	VPG-BTD	\$55
Go Direct Projectile Launcher	GDX-PL	\$499
Vernier Projectile Launcher	VPL	\$440
Projectile Stop	PS-VPL	\$40
Pulley Bracket	B-SPA	\$12
Go Direct Rotary Motion	GDX-RMS ⚡	\$195
Rotary Motion Sensor	RMV-BTD ⚡	\$189
Rotational Motion Accessory Kit	AK-RMV ⚡	\$129
Time of Flight Pad	TOF-VPL	\$84
Ultra Pulley Attachment	SPA	\$24

## Waves and Sound

Product	Order Code	Price
Microphone	MCA-BTA	\$55
Sound Level Sensor	SLS-BTA	\$69

## Thermodynamics

Product	Order Code	Price
Gas Pressure Sensor	GPS-BTA	\$103
Stainless Steel Temperature Probe	TMP-BTA	\$45
Surface Temperature Sensor	STS-BTA	\$34

## Electricity and Magnetism

Product	Order Code	Price
Magnetic Field Sensor	MG-BTA	\$69
Power Amplifier	PAMP	\$249
Differential Voltage Probe	DVP-BTA	\$49
Current Probe	DCP-BTA	\$49
Instrumentation Amplifier	INA-BTA	\$89
Optional Breadboard Kit for the Vernier Circuit Board 2	VCB2-0BBK	\$29
Extech® Digital DC Power Supply	EXPS	\$259
Charge Sensor	CRG-BTA	\$89
Solenoid	PHY-SN700	\$79

## Light and Optics

Product	Order Code	Price
Light Sensor	LS-BTA ⚡	\$69
Polarizer/Analyzer Set for Optics Expansion Kit	PAK-0EK	\$89
Combination 1.2 m Track/Optics Bench	TRACK	\$160
Combination 2.2 m Track/Optics Bench	TRACK-LONG	\$272
Green Diffraction Laser	GDL-DAK	\$250

## Lab Books

Product	Order Code	Price
<i>Physics with Vernier</i>	PWV*	\$59
<i>Advanced Physics with Vernier—Mechanics</i>	PHYS-AM*	\$59
<i>Advanced Physics with Vernier—Beyond Mechanics</i>	PHYS-ABM*	\$59
<i>Physics Explorations and Projects</i>	PEP*	\$59
<i>Vernier Video Analysis: Motion and Sports</i>	HSB-VVAMS-E	\$29
<i>Vernier Video Analysis: Conservation Laws and Forces</i>	HSB-VVACLF-E	\$29

\* Includes printed book and download; also available as a download only

To see the full suite of Vernier physics products, please visit [vernier.com/physics](https://www.vernier.com/physics)



# Engineering

## Why Vernier?

Vernier engineering solutions harness the power of analytical software and the precision of high-quality sensors to help students sharpen their design skills and prepare to enter the workforce. As with all of our solutions, our engineering technology is backed by unparalleled support.

### Quality

Durable hardware for lab and field use

### Affordable

Designed for education and education budgets

### Versatile

Supports a variety of devices and experiments



“ Our projects are about more than just supporting the need for engineering education in local classrooms. Vernier products help deepen our students' learning through experiential, hands-on community engagement.

*Maija A. Benitz, PhD  
Associate Professor of Engineering  
Roger Williams University*

[vernier.com/engineering](https://www.vernier.com/engineering)

# A Guide to Vernier Data Collection

We're here to support you as an educator as you incorporate data-collection technology into your instruction. See how our products provide you with affordable laboratory solutions designed for student success.

Our Guarantee: Most of our products are protected by a 5-year limited warranty. And after 5 years? We'll make every attempt to repair your equipment.

## What You Need to Get Started with Go Direct Sensors

Go Direct® sensors connect directly (no interface required) to your computer, Chromebook™, or compatible mobile device via USB or Bluetooth® wireless technology.

- A Go Direct sensor**
- B Computer, Chromebook, tablet, smartphone, LabQuest 3**
- C Software**
  - Vernier Graphical Analysis® Pro
  - Arduino IDE
  - National Instruments LabVIEW™
  - Python®
  - JavaScript™

## What You Need to Get Started with LabQuest Sensors

LabQuest® sensors have a cable with a plug that makes it easy to connect to an interface without any additional wiring. Use LabQuest sensors with a Vernier interface, Arduino®, NI ELVIS, NI myDAQ, or your own DAQ hardware. If using a non-Vernier interface, these sensors require a +5.0 volt supply voltage and output a 0 to 5 volt signal. Most sensors have a simple, linear calibration.

- A LabQuest sensor**
- B Interface (LabQuest, DAQ, or Arduino)**
- C Computer, Chromebook, tablet, smartphone**
- D Software**
  - Vernier Graphical Analysis Pro
  - National Instruments LabVIEW
  - Arduino IDE
  - Python

# Vernier Graphical Analysis Pro

Our award-winning app enables advanced data analysis, provides opportunities to reinforce and extend learning, and includes features that support hybrid teaching models. We're always working to add and improve features—check out the “What’s New” tab in the app.


**1-Year Site License (unlimited seats)** GAP-1YR \$149

**3-Year Site License (unlimited seats)** GAP-3YR \$299

### Free Trial for Educators

Try out Graphical Analysis Pro for free for 30 days. Access the sample experiments and enhanced analysis tools to use with your students.

Get a free trial and learn about site license options at [vernier.com/graphical-analysis](https://www.vernier.com/graphical-analysis)

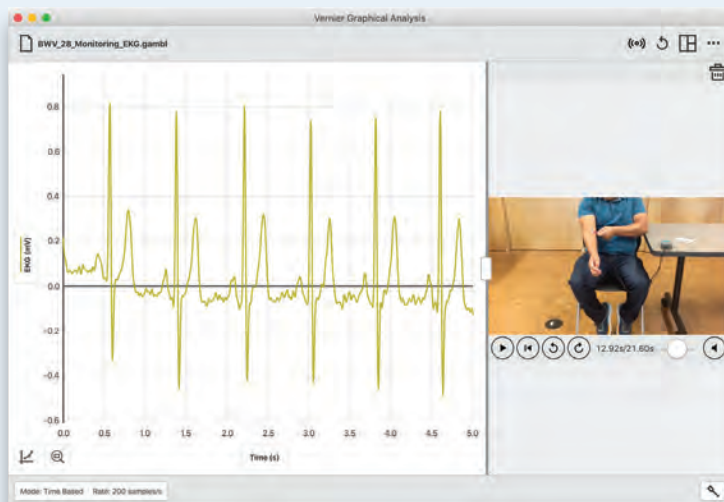

Get everything you need with our new 4-in-1 software bundle!  
**3-Year Site License** DCA-3YR \$499 | Learn more at [vernier.com/dca-3yr](https://www.vernier.com/dca-3yr)

## Free 30-Day Preview

Not sure if our technology is right for your laboratory? Our risk-free preview program makes it easier to decide. Most Vernier products are available for a 30-day preview (or longer, if requested) to educational institutions.


# Introduction to Engineering

## Go Direct Structures & Materials Tester

Use our Go Direct Structures & Materials Tester to evaluate the strength of model bridges and engineered structures by measuring the applied load. Utilizing both load and displacement sensors, your students can evaluate the properties of materials.

GDX-VSMT \$1,199

[vernier.com/gdx-vsmt](http://vernier.com/gdx-vsmt)

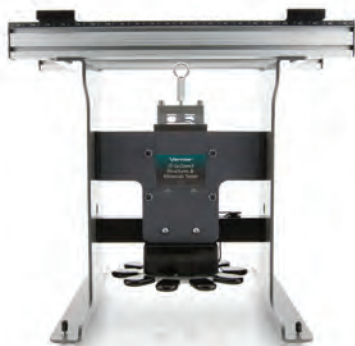

## Advanced Wind Experiment Kit

Use this kit to teach the engineering design process. Investigate different blade designs, gear ratios, and generators.

KW-AWX \$174

[vernier.com/kw-awx](http://vernier.com/kw-awx)



## Vernier Variable Load

Use the Vernier Variable Load in conjunction with Go Direct Energy to provide a range of resistive loads for projects such as engineering wind turbines or investigating solar panels. Students can adjust the potentiometer to provide resistances between 2 and 200  $\Omega$  to determine the optimal load on a system.

VES-VL \$79 [vernier.com/ves-vl](http://vernier.com/ves-vl)



## Go Direct Energy

Go Direct Energy measures voltage and current as well as displays power and energy output of scale model wind turbines and solar panels, so students can quantitatively evaluate the effects of their design changes.

GDX-NRG \$99 [vernier.com/gdx-nrg](http://vernier.com/gdx-nrg)



See all our products for engineering at [vernier.com/engineering](http://vernier.com/engineering)

# Wind Turbine Design

## Tower and Base Set

Do you need a tower for your turbine nacelle? This is the same tower that comes in the Advanced Wind Experiment Kit. The tower has a diameter that fits inside 1-inch PVC fittings.

KW-TBS \$29

[vernier.com/kw-tbs](http://vernier.com/kw-tbs)



## Basic Turbine Building Parts

The Basic Turbine Building Parts kit includes three hubs, a wind turbine generator, and 25 dowels, all in one package.

KW-BTPART \$16

[vernier.com/kw-btpart](http://vernier.com/kw-btpart)



## Wind Turbine Generator with Wires

This is the primary generator for wind turbine experiments because it runs smoothly and provides high power output at a relatively low RPM.

KW-GEN3 \$20

[vernier.com/kw-gen](http://vernier.com/kw-gen)



## Wind Turbine Hub (3 Pack)

With these 12-hole crimping hubs, made from recycled plastic, students can turn a DC generator into a wind turbine.

KW-WTH3 \$22

[vernier.com/kw-wth3](http://vernier.com/kw-wth3)



## Nacelle

Build your own wind turbine and use the KidWind Nacelle to hold the generator (motor). Package includes the KidWind Nacelle and the Motor Mount Pack. You will also need a generator and a way to affix the turbine blades.

KW-NAC \$18

[vernier.com/kw-nac](http://vernier.com/kw-nac)



## simpleGEN

Students can use the easy-to-build AC generator of the simpleGEN to explore the basics of electrical generator design.

KW-SGEN \$79

[vernier.com/kw-sgen](http://vernier.com/kw-sgen)



## Gear Set

The small 8-tooth gear fits on 2 mm driveshafts that are found on many DC generators. The gears have a keying feature and can be changed quickly and easily using the included hex locks. The hex locks secure to our 9/64 inch hex driveshaft, which is included in the Drivetrain Set (KW-DS, \$8).

Gear sizes: 64 teeth, 32 teeth, 16 teeth, 8 teeth

KW-GEAR \$9

[vernier.com/kw-gear](http://vernier.com/kw-gear)



## Photogate Bracket

The Photogate Bracket enables the measurement of rotations per minute (RPM) of a wind turbine with a Vernier photogate.

KW-PGBRAC \$12

[vernier.com/kw-pgbrac](http://vernier.com/kw-pgbrac)





## Biomedical Engineering with Go Direct Sensors

With wireless options and multiple on-board sensors, Go Direct® sensors are perfect for analyzing and studying physiological functions.

### Go Direct EKG

Go Direct EKG has five channels: EKG, heart rate, EMG, EMG rectified, and voltage.

GDX-EKG \$175

[vernier.com/gdx-ekg](http://vernier.com/gdx-ekg)



### Go Direct O<sub>2</sub> Gas

This sensor measures gaseous oxygen concentration levels and air temperature.

GDX-O2 \$219

[vernier.com/gdx-o2](http://vernier.com/gdx-o2)



### Go Direct Blood Pressure

Go Direct Blood Pressure has seven channels: cuff pressure, mean arterial pressure, systolic pressure, diastolic pressure, pulse rate, oscillations, and envelope.

GDX-BP \$125

[vernier.com/gdx-bp](http://vernier.com/gdx-bp)



### Go Direct Force Plate

Go Direct Force Plate has six channels: force, hang time, jump height, X balance, Y balance, and auxiliary voltage.

GDX-FP \$399

[vernier.com/gdx-fp](http://vernier.com/gdx-fp)



### Go Direct Spirometer

Go Direct Spirometer has six channels: flow rate, volume, adjusted volume, cycle volume, respiration rate, and differential pressure.

GDX-SPR \$229

[vernier.com/gdx-spr](http://vernier.com/gdx-spr)



### Go Direct Surface Temperature

With a range of -25 to 125°C, this sensor is designed for use in situations in which low thermal mass or flexibility is required, such as on human skin.

GDX-ST \$99

[vernier.com/gdx-st](http://vernier.com/gdx-st)



### Go Direct Hand Dynamometer

Go Direct Hand Dynamometer has seven channels: force, x-axis acceleration, y-axis acceleration, z-axis acceleration, x-axis gyro, y-axis gyro, and z-axis gyro.

GDX-HD \$135

[vernier.com/gdx-hd](http://vernier.com/gdx-hd)



### Go Direct Respiration Belt

Go Direct Respiration Belt has four channels: force, respiration rate, steps, and step rate.

GDX-RB \$125

[vernier.com/gdx-rb](http://vernier.com/gdx-rb)



### Go Direct CO<sub>2</sub> Gas

Go Direct CO<sub>2</sub> Gas has three channels: CO<sub>2</sub> gas, temperature, and relative humidity.

GDX-CO2 \$229

[vernier.com/gdx-co2](http://vernier.com/gdx-co2)



### Go Direct Acceleration

This 3-axis acceleration sensor has two acceleration ranges ( $\pm 157 \text{ m/s}^2$  and  $\pm 1960 \text{ m/s}^2$ ) plus an altimeter and a 3-axis gyroscope.

GDX-ACC \$115

[vernier.com/gdx-acc](http://vernier.com/gdx-acc)



# Arduino with Vernier Sensors


Help students build coding skills and foster creative critical thinking by combining Vernier sensors with Arduino.®

## Vernier Coding Activities with Arduino®: Analog Sensors

Combining Vernier sensors with easy-to-program Arduino microcontrollers is an engaging way for students to make connections between coding and the physical world. The 8 lessons in this e-book introduce students to programming functions, logic statements, and more.

DOWNLOAD ONLY  
VCA-AS-E \$24

[vernier.com/vca-as-e](http://vernier.com/vca-as-e)

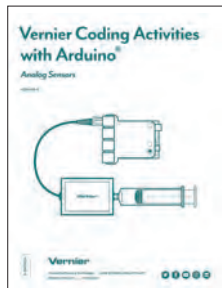

## Vernier Arduino Interface Shield

Conveniently connect the SparkFun® RedBoard or Arduino Uno to Vernier LabQuest® sensors with the Vernier Arduino Interface Shield.

BT-ARD \$29

[vernier.com/bt-ard](http://vernier.com/bt-ard)



## Gas Pressure Sensor

Use the Gas Pressure Sensor to monitor pressure changes of a gas.

GPS-BTA \$103 [vernier.com/gps-bta](http://vernier.com/gps-bta)



## Surface Temperature Sensor

Measure temperature where low thermal mass or flexibility is required.

Range: -25 to 125°C

STS-BTA \$34 [vernier.com/sts-bta](http://vernier.com/sts-bta)



## Motion Detector

The Motion Detector uses ultrasound to measure the position of objects.

Range: 0.15 to 6 m

Resolution: 1 mm

MD-BTD \$105 [vernier.com/md-btd](http://vernier.com/md-btd)



## pH Sensor

This is a general-purpose pH sensor.

Range: pH 0 to 14

Accuracy: ±0.2 pH units

PH-BTA \$105 [vernier.com/ph-bta](http://vernier.com/ph-bta)



## Digital Control Unit

Use the digital output lines of an interface to control DC electrical devices.

DCU-BTD \$69 [vernier.com/dcu-btd](http://vernier.com/dcu-btd)



## Anemometer

This is an impeller-type anemometer for measuring wind speed.

Range: 0.5 to 30 m/s (1 to 67 mph)

ANM-BTA \$109 [vernier.com/anm-bta](http://vernier.com/anm-bta)



Read the online guide and see all our products for Arduino at [vernier.com/arduino](http://vernier.com/arduino)

## Sensors and Coding

Students gain a deeper understanding of coding when they see their code interact with sensor data. Our coding solutions give students valuable experience using popular programming languages to affect things beyond the screen.



### LabVIEW

Use our free activity book to introduce NI LabVIEW™ programming with Vernier sensors in your classroom.

[vernier.com/ni-labview](http://vernier.com/ni-labview)

# National Instruments LabVIEW and Vernier

Introduce your students to NI LabVIEW software, a programming language used throughout the engineering disciplines. We have sample LabVIEW programs (VIs) for LabQuest Mini, myDAQ, Go Direct® sensors, and other Vernier hardware.

## With LabQuest Sensors



### myDAQ Adapter

The myDAQ Adapter can be used to perform data acquisition with more than 75 Vernier LabQuest sensors and the NI myDAQ interface (sold separately). It is designed for use with NI LabVIEW software.



BT-MDAQ \$69

[vernier.com/bt-mdaq](http://vernier.com/bt-mdaq)

### Analog Protoboard Adapter

Use this adapter to connect Vernier LabQuest sensors to a non-Vernier interface, such as NI ELVIS. The connector fits into a standard prototyping board.



BTA-ELV \$12

[vernier.com/bta-elv](http://vernier.com/bta-elv)

### LabQuest Mini

LabQuest Mini is a powerful, affordable, and easy-to-use sensor interface for data acquisition with more than 75 Vernier LabQuest sensors.

LQ-MINI \$189

[vernier.com/lq-mini](http://vernier.com/lq-mini)

## With Go Direct Sensors

Integrate over 50 wireless sensors into your LabVIEW project to acquire data or control your NI DAQ hardware.

### Go Direct Acceleration



GDX-ACC \$115 [vernier.com/gdx-acc](http://vernier.com/gdx-acc)

### Go Direct Motion



GDX-MD \$129 [vernier.com/gdx-md](http://vernier.com/gdx-md)

### Go Direct Force and Acceleration



GDX-FOR \$119 [vernier.com/gdx-for](http://vernier.com/gdx-for)

### Go Direct Light and Color



GDX-LC \$99 [vernier.com/gdx-lc](http://vernier.com/gdx-lc)

### Go Direct Rotary Motion



GDX-RMS ⚡ \$195 [vernier.com/gdx-rms](http://vernier.com/gdx-rms)

### Go Direct Weather



GDX-WTHR \$135 [vernier.com/gdx-wthr](http://vernier.com/gdx-wthr)

See all our products for NI LabVIEW at [vernier.com/ni-labview](http://vernier.com/ni-labview)


⚡ For California Proposition 65 warning, see page 36.



## Python and VPython

We have created a guide and examples that make it easy to start using Vernier sensors with Python® and VPython.

[vernier.com/python](http://vernier.com/python)

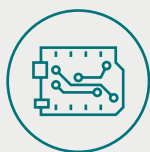

## Arduino IDE

Use our starter package to teach the basics of coding Arduino with Vernier sensors.

[vernier.com/arduino](http://vernier.com/arduino)



## JavaScript

Students can use JavaScript™ to integrate Go Direct sensor data into their custom web applications.

[vernier.com/javascript](http://vernier.com/javascript)

# Engineering Products

## Wind Turbine Design

Product	Order Code	Price
Advanced Wind Experiment Kit	KW-AWX	\$174
Balsa Blade Sheets (10 Sheets)	KW-BBS10	\$12
Basic Turbine Building Parts	KW-BTPART	\$16
Drivetrain Set	KW-DS	\$8
Gear Set	KW-GEAR	\$9
Go Direct® Energy	GDX-NRG	\$99
Hub (3 Pack)	KW-WTH3	\$22
KidWind Photogate Bracket	KW-PGBRAC	\$12
Nacelle	KW-NAC	\$18
Tower and Base Set	KW-TBS	\$29
Vernier Variable Load	VES-VL	\$79
Wind Turbine Generator with Wires	KW-GEN3	\$20

## Engineering with Arduino

Product	Order Code	Price
Anemometer	ANM-BTA	\$109
Digital Control Unit	DCU-BTD	\$69
Gas Pressure Sensor	GPS-BTA	\$103
Motion Detector	MD-BTD	\$105
pH Sensor	PH-BTA	\$105
Surface Temperature Sensor	STS-BTA	\$34
Vernier Arduino Interface Shield	BT-ARD	\$29
<i>Vernier Coding Activities with Arduino:</i> <i>Analog Sensors lab book</i>	VCA-AS-E	\$24

## LabQuest Sensors

Sensor	Order Code	Price
Barometer	BAR-BTA	\$79
Gas Pressure Sensor	GPS-BTA	\$103
Light Sensor	LS-BTA ⚡	\$69
Magnetic Field Sensor	MG-BTA	\$69
Microphone	MCA-BTA	\$55
Soil Moisture Sensor	SMS-BTA	\$148
Stainless Steel Temperature Probe	TMP-BTA	\$45

Learn more about over 80 LabQuest sensors at [vernier.com/labquest](https://www.vernier.com/labquest)

## Biomedical Engineering

Sensor	Order Code	Price
Go Direct Acceleration	GDX-ACC	\$115
Go Direct Blood Pressure	GDX-BP	\$125
Go Direct CO <sub>2</sub> Gas	GDX-CO2	\$229
Go Direct EKG	GDX-EKG	\$175
Go Direct Force Plate	GDX-FP	\$399
Go Direct Hand Dynamometer	GDX-HD	\$135
Go Direct O <sub>2</sub> Gas	GDX-O2	\$219
Go Direct Respiration Belt	GDX-RB	\$125
Go Direct Spirometer	GDX-SPR	\$229
Go Direct Surface Temperature	GDX-ST	\$99

## NI LabVIEW and Vernier

Product	Order Code	Price
Analog Protoboard Adapter	BTA-ELV	\$12
myDAQ Adapter	BT-MDAQ	\$69
LabQuest® Mini	LQ-MINI	\$189

## Go Direct Sensors

Sensor	Order Code	Price
Go Direct Acceleration	GDX-ACC	\$115
Go Direct Force and Acceleration	GDX-FOR	\$119
Go Direct Light and Color	GDX-LC	\$99
Go Direct Motion	GDX-MD	\$129
Go Direct Rotary Motion	GDX-RMS ⚡	\$195

Learn more about over 65 Go Direct sensors at [vernier.com/go-direct](https://www.vernier.com/go-direct)

See all of our engineering products online at [vernier.com/engineering](https://www.vernier.com/engineering)



# Index

## A

ADInstruments 2  
Advanced Chemistry with Vernier 13  
Advanced Wind Experiment Kit 29  
Analog Protoboard Adapter 33  
Anemometer 32  
Arduino® IDE 33  
Arduino Interface Shield 32

## B

Balances [vernier.com/ohaus](http://vernier.com/ohaus)  
BioChamber 250 [vernier.com/bc-250](http://vernier.com/bc-250)  
BioChamber 2000 [vernier.com/bc-2000](http://vernier.com/bc-2000)  
Biology Go Direct Standard Package 5  
Biology with Vernier 5  
Bio-Rad® 7  
BlueView Transilluminator 7

## C

Cart Fan 21  
Celestron® Digital Microscope Imagers  
[vernier.com/cs-5mp](http://vernier.com/cs-5mp)  
Charging stations  
Go Direct [vernier.com/gdx-crg](http://vernier.com/gdx-crg)  
Go Direct Sensor Cart [vernier.com/gdx-cart-crg](http://vernier.com/gdx-cart-crg)  
LabQuest 3 [vernier.com/lq3-crg](http://vernier.com/lq3-crg)  
Circuit Board 23  
Color Mixer Kit 24

## D

Diffraction Apparatus 24  
Digital Control Unit 32  
Dynamics Cart and Track System with  
Go Direct Sensor Cart 21

## E

Electrochemistry Half-Cell Plate 16  
Electrochemistry Metals Kit 16  
Electrode amplifiers  
Electrode Amplifier [vernier.com/ea-bta](http://vernier.com/ea-bta)  
Go Direct Electrode Amplifier  
[vernier.com/gdx-ea](http://vernier.com/gdx-ea)  
Ion-Selective Electrode Amplifier  
[vernier.com/gdx-isea](http://vernier.com/gdx-isea)  
Electrostatic High-Voltage Genecon 23  
Electrostatics kits  
Electrostatic Kit 23  
High-Voltage Electrostatics Kit 23  
ELVIS protoboard adapters  
[vernier.com/protoboard-adapters](http://vernier.com/protoboard-adapters)  
Emissions spectrometer 25  
Environmental Science LabQuest 3 Package 9  
Equipment return 37  
Extech® Digital Power Supply 23

## F

Frequency Generator 22

## G

Gas Chromatograph 15  
Gas pressure sensors  
Gas Pressure Sensor 32  
Go Direct Gas Pressure 13, 22  
Go Direct Wide-Range Pressure 16  
Pressure Sensor 400  
[vernier.com/ps400-bta](http://vernier.com/ps400-bta)  
Go Direct Centripetal Force System 21  
Go Direct Sensor Clamp 9

Go Direct sensors

Go Direct 3-Axis Magnetic Field 23  
Go Direct Acceleration 31, 33  
Go Direct Blood Pressure 6, 31  
Go Direct CO<sub>2</sub> Gas 5, 31  
Go Direct Conductivity 7, 8  
Go Direct Current 23  
Go Direct Cyclic Voltammetry System 16  
Go Direct Drop Counter 13  
Go Direct EKG 6, 31  
Go Direct Energy 29  
Go Direct Force and Acceleration 21, 33  
Go Direct Force Plate 21, 31  
Go Direct Gas Pressure 13, 22  
Go Direct Hand Dynamometer 6, 31  
Go Direct Light and Color 24, 33  
Go Direct Melt Station 15  
Go Direct Mini GC 15  
Go Direct Motion 21, 33  
Go Direct Nitrate Ion-Selective Electrode 8  
Go Direct O<sub>2</sub> Gas 6, 31  
Go Direct Optical Dissolved Oxygen 5, 8  
Go Direct ORP [vernier.com/gdx-orp](http://vernier.com/gdx-orp)  
Go Direct PAR 8  
Go Direct pH 13  
Go Direct Photogate 21  
Go Direct Platinum-Cell Conductivity  
[vernier.com/gdx-compt](http://vernier.com/gdx-compt)  
Go Direct Polarimeter 14  
Go Direct Pyranometer [vernier.com/gdx-pyr](http://vernier.com/gdx-pyr)  
Go Direct Radiation Monitor 25  
Go Direct Respiration Belt 6, 31  
Go Direct Rotary Motion 33  
Go Direct Salinity 8  
Go Direct Soil Moisture 8  
Go Direct Sound 22  
Go Direct Spectrometers 5, 7, 13–17, 25  
Go Direct Spirometer 6, 31  
Go Direct Static Charge 23  
Go Direct Structures & Materials Tester 29  
Go Direct Surface Temperature 22, 31  
Go Direct Temperature 13, 22  
Go Direct Thermocouple [vernier.com/gdx-tc](http://vernier.com/gdx-tc)  
Go Direct Tris-Compatible Flat pH 5, 7, 8, 14  
Go Direct Voltage 23  
Go Direct Weather Sensor 33  
Go Direct Weather System 8  
Go Direct Wide-Range Pressure 17  
Go Direct Wide-Range Temperature 15

Graphical Analysis Pro app 4, 12, 20, 28

Green Diffraction Laser 24

## H

Human Physiology Experiments: Volume 1 6  
Human Physiology Experiments: Volume 2 6  
Human Physiology Go Direct Standard Package 6

## I

Instrumental Analysis app 15  
International sales 37  
Investigating Biology through Inquiry 5  
Investigating Environmental Science  
through Inquiry 9

## J

JavaScript™ 33

## L

LabArchives 2  
LabQuest 3 4, 12, 20  
LabQuest Mini 33  
LabQuest Viewer [vernier.com/lq-view](http://vernier.com/lq-view)  
LabVIEW™ 32

## M

Mirror Set 24  
Moment of Inertia Accessory Kit 21  
Motion Detector 32  
myDAQ Adapter 33

## O

Optics Expansion Kit 24  
Organic Chemistry with Vernier 15

## P

pH Sensor 32  
pH Storage Solution [vernier.com/ph-ss](http://vernier.com/ph-ss)  
Polarimeter 14  
Polarizer/Analyzer Set 24  
Power Amplifier 22  
Power Amplifier Accessory Speaker 22  
Precision Volume Dispenser 16  
Primary Productivity Kit [vernier.com/ppk](http://vernier.com/ppk)  
Privacy policy 37  
Prop 65 (California) 36  
Python® and VPython 33

## Q

Qubit Systems sensors [vernier.com/qubit](http://vernier.com/qubit)

## R

Radiation Monitor 25  
Renewable Energy with Vernier 9  
Resonance Apparatus 22  
Returns 37

## S

Software licenses 37  
Spectral Analysis app 7, 14, 25  
Spectrometers/Spectrophotometers  
Go Direct Emissions Spectrometer 25  
Go Direct Fluorescence/UV-VIS  
Spectrophotometer 7, 14, 17  
Go Direct SpectroVis Plus 5, 7, 13  
Go Direct UV-VIS Spectrophotometer 7, 15  
Go Direct Visible Spectrophotometer 17  
Vernier Emissions Fiber 25  
Vernier Flash Photolysis Spectrometer 17  
Vernier Spectrophotometer Optical Fiber 16  
Spectrum Tube Power Supply 25  
Spectrum Tubes 25  
Stir Station 13  
Surface Temperature Sensor 32

## T

Technical support 37  
Texas Instruments software  
[vernier.com/ti-software](http://vernier.com/ti-software)

## V

Variable Load 29  
Vernier Coding Activities with Arduino®:  
Analog Sensors 32  
Video Analysis app 20  
Vernier Circuit Board 2 23

## W

Warranty information 37  
Water Depth Sampler [vernier.com/wds](http://vernier.com/wds)  
Water Quality Bottles [vernier.com/wq-bot](http://vernier.com/wq-bot)  
Water Quality with Vernier 9  
Wind turbine design products 30

# California Proposition 65 Warning

📌 **PROP 65**—For more information, go to [P65Warnings.ca.gov](https://www.vernier.com/p65warnings)

Vernier Products Affected	WARNING
BioChamber 2000	⚠️ WARNING: This product can expose you to chemicals, including bisphenol A (BPA), which are known to the State of California to cause cancer, and methyl isobutyl ketone (MIBK), which is known to the State of California to cause birth defects or other reproductive harm.
BioChamber 250	⚠️ WARNING: This product can expose you to chemicals, including methyl isobutyl ketone (MIBK), which are known to the State of California to cause cancer and birth defects or other reproductive harm.
BlueView Transilluminator	⚠️ WARNING: This product can expose you to chemicals, including ethyl acrylate, which are known to the State of California to cause cancer.
Electrochemistry Metals Kit	⚠️ WARNING: This product can expose you to chemicals, including lead, which are known to the State of California to cause cancer.
Electrostatics Kit	⚠️ WARNING: This product can expose you to chemicals, including nickel, which are known to the State of California to cause cancer.
Go Direct® Melt Station	⚠️ WARNING: This product can expose you to chemicals, including nickel, which are known to the State of California to cause cancer.
Go Direct Rotary Motion Sensor	⚠️ WARNING: This product can expose you to chemicals, including chromium, which are known to the State of California to cause cancer and birth defects or other reproductive harm.
Hanging Mass Set 250 g	⚠️ WARNING: This product can expose you to chemicals, including lead, which are known to the State of California to cause cancer.
Human Physiology Go Direct Standard Package	⚠️ WARNING: This product can expose you to chemicals, including methyl isobutyl ketone (MIBK), which are known to the State of California to cause cancer and birth defects or other reproductive harm.
Light Sensor	⚠️ WARNING: This product can expose you to chemicals, including antimony, which are known to the State of California to cause cancer.
Melt Station	⚠️ WARNING: This product can expose you to chemicals, including nickel, which are known to the State of California to cause cancer.
OHAUS Scout® 120 g	⚠️ WARNING: This product can expose you to chemicals, including Di(2-ethylhexyl) phthalate (DEHP), which are known to the State of California to cause cancer and birth defects or other reproductive harm.
OHAUS Scout 220 g	⚠️ WARNING: This product can expose you to chemicals, including Di(2-ethylhexyl) phthalate (DEHP), which are known to the State of California to cause cancer and birth defects or other reproductive harm.
OHAUS Scout 420 g	⚠️ WARNING: This product can expose you to chemicals, including Di(2-ethylhexyl) phthalate (DEHP), which are known to the State of California to cause cancer and birth defects or other reproductive harm.
Power Amp Accessory Speaker	⚠️ WARNING: This product can expose you to chemicals, including chromium, which are known to the State of California to cause cancer and birth defects or other reproductive harm.
Rotary Motion Motor Kit	⚠️ WARNING: This product can expose you to chemicals, including formaldehyde, which are known to the State of California to cause cancer.
Rotational Motion Accessory Kit	⚠️ WARNING: This product can expose you to chemicals, including chromium, which are known to the State of California to cause cancer and birth defects or other reproductive harm.
Spectrum Tube Power Supply	⚠️ WARNING: This product can expose you to chemicals, including nickel, which are known to the State of California to cause cancer.
Spectrum Tubes (Air, Argon, Carbon Dioxide, Hydrogen, Helium, Neon, Nitrogen)	⚠️ WARNING: This product can expose you to chemicals, including nickel, which are known to the State of California to cause cancer.
Vernier Circuit Board 2	⚠️ WARNING: This product can expose you to chemicals, including nickel, which are known to the State of California to cause cancer.
Vernier Rotary Motion Sensor	⚠️ WARNING: This product can expose you to chemicals, including chromium, which are known to the State of California to cause cancer and birth defects or other reproductive harm.

# Doing Business with Us

## Satisfaction Guarantee

Vernier has been selling science education software and data-collection hardware since 1981. We pride ourselves on the quality and affordability of our products and our service to our customers. If at any time you are unhappy with any of our products or service, please get in touch.

### Vernier Science Education

13979 SW Millikan Way  
Beaverton, OR 97005-2886  
vernier.com • info@vernier.com  
Toll Free: 888-VERNIER (888-837-6437)  
Fax: 503-277-2440

## Product Usage

Vernier products are designed for educational use. Our products are not designed nor are they recommended for any industrial, medical, or commercial process, such as life support, patient diagnosis, control of a manufacturing process, or industrial testing of any kind. We design our products with the specifications and features that educators and students need to be successful. In our effort to keep our products affordable and easy to use, we may not meet the specifications or include the features that an industrial scientist or medical professional might want.

## Equipment Return

Any product that does not meet your needs may be returned within 30 days for a full refund. Equipment returned after 30 days may be subject to a restocking fee.

A Return Merchandise Authorization, available from Vernier, is required for any product return. Equipment returned for exchange or credit must be in new condition and in its original packaging.

## Prices and Shipping

Prices are effective January 1, 2025 and supersede previously published prices. Prices are in US dollars and are FOB shipping point. Shipping charges may vary, depending on method of shipment. Increased shipping charges for heavier or bulkier items may apply due to weight or dimensions. Applicable sales tax may be charged. Prices in this catalog are for US educational institutions only and are subject to change without notice.

## Preview Policy

Most Vernier products are available for a 30-day preview (or longer, if requested) to US educational institutions.

## International Sales

All Vernier orders for use outside of the US and Canada are handled by us and the worldwide network of Vernier dealers. Contact us for more information.

Sales of Vernier products in Canada are handled by

### Vernier Canada

7030 Woodbine Ave. Suite 500  
Markham, Ontario L3R 6G2  
Canada  
verniercanada.ca • info@verniercanada.ca  
Phone: (800) 376-4210 • Local: (705) 915-3656

## Warranties

Most Vernier-branded products carry a 5-year limited warranty. Product-specific details can be found under the Support tab on each product's web page. During the warranty period, Vernier will repair or replace the item if there is a defect in materials or workmanship. Outside the warranty, Vernier will attempt to repair most products. The Vernier warranty covers products when used by educational institutions only. Products manufactured by anyone other than Vernier are subject to the conditions of the warranty supplied by the manufacturer.

Additional exclusions and limitations can be found at [vernier.com/warranty](http://vernier.com/warranty)

## Software Licenses

Vernier Graphical Analysis, Vernier Spectral Analysis, and Vernier Instrumental Analysis are available as free downloads from our website, distributed as a progressive web app, or distributed through the appropriate web store. Vernier Graphical Analysis Pro is available as a subscription service. Vernier Video Analysis is available as a subscription service and is distributed as a progressive web app. Apps for iOS, iPadOS, Android, and Chrome are distributed through their respective stores. Terms and licensing are thus determined entirely by these stores.

## Other Software

Software from other companies, such as Texas Instruments, is licensed under separate agreements by their respective companies.

## Privacy Policy

Vernier Science Education does not sell, lease, or loan our mailing list or portions thereof to anyone at any time. We do not store credit card information on our online store or in our accounting system. For more information on our privacy policy, see [vernier.com/legal](http://vernier.com/legal)

If you wish to be removed from our mailing list, simply write to us at [updates@vernier.com](mailto:updates@vernier.com), and we will remove you immediately.

## Trademarks

Vernier Science Education, LabQuest, SpectroVis, Vernier and caliper design, Go Direct, Go Wireless, LabQuest Viewer, Vernier Graphical Analysis, Vernier Spectral Analysis, Vernier Video Analysis, and Vernier Instrumental Analysis are our registered trademarks. Vernier.com and BlueView are our trademarks or trade dress.

Apple, the Apple logo, iPhone, iPad, iPadOS, and macOS are trademarks of Apple Inc., registered in the United States and other countries. App Store is a service mark of Apple Inc.

Arduino® and  are trademarks of Arduino SA.

National Instruments, NI, and LabVIEW are trademarks or trade names of National Instruments Corporation.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Vernier Science Education is under license.

All other marks not owned by us that appear herein are the property of their respective owners, who may or may not be affiliated with, connected to, or sponsored by us.

## Technical Support

We are readily available to help you with individual questions about our software and hardware—simply email [support@vernier.com](mailto:support@vernier.com), chat with us live on [vernier.com](http://vernier.com), or call us at our toll free number: 888-VERNIER (888-837-6437).

Our email newsletter makes it easy to access new ideas, learn about new products, and get inspired by fellow educators. Sign up at [vernier.com/newsletter](http://vernier.com/newsletter)

## Legal

Visit our Legal Center at [vernier.com/legal](http://vernier.com/legal) to find our privacy statements, terms of use, and other information about our products and services.



Vernier Science Education  
13979 SW Millikan Way  
Beaverton, OR 97005-2886

888-VERNIER  
(888-837-6437)  
fax 503-277-2440

vernier.com  
info@vernier.com



Recipient not at your school?  
Please send updates to  
updates@vernier.com

## Try Vernier Products for Free

Not sure if our technology is right for your lab? Our risk-free preview program makes it easy to decide. Most Vernier products are available for a 30-day preview (or longer, if requested) to educational institutions.

## Our Guarantee

Most of our products are protected by a 5-year limited warranty. And after 5 years? We'll make every attempt to repair your equipment.